



Shenchen Pump

We aim to build the world brand of Peristaltic Pump

Exported to **63** countries



Peristaltic Pump • Syringe Pump

Baoding Shenchen Precision Pump Co.,Ltd

Baoding Peristaltic Pump Engineering Technology Research Center

(2016.2)

SHENCHEN

Built in 2006 (2006-2008)

- Innovation—the first producer of easy load pump head with PPS material, solving the technical problem what the PSU material pump head can't bearing corrosion resistance.
 - Successfully launch the DG series pump head
 - Successfully launch the M series standard peristaltic pump, K series flow rate peristaltic pump, FJ series dispensing peristaltic pump.
- Shenchen products recognized by the market quickly with rich product series and excellent quality.

2006

In 2009

- Successfully launch MC series pump head, solve the technical problems of inconvenient card of traditional DG series multi channel pump head, the design is classic work of micro innovation.
 - Successfully launch MiniPump series compact pump head, SN series standard pump head.
 - Successfully launch N series standard peristaltic pump.
- Shenchen products continuously fill gaps in the market.

2009

2010

In 2010-2013

- Successfully registered "shenchen" trademark
- Successfully got the CE certification from European Union
- Successfully launch the intelligent filling system CF600, DF600. The filling accuracy is high to 0.5%, obtain the software copyright, create a new era of peristaltic pump filling application.
- Successfully launch the V series, F series intelligent type peristaltic pump, bring customers new using experience with the high-quality color animation interface, touch sensitive effect, and high accuracy.
- Successfully launch SPLab series syringe pump
- Shenchen move to the provincial key industrial park.
- Shenchen upgrade all the products and successfully entered the international market.

2016

In 2014-2016

- Shenchen launch the LabN, LabV, LabF series of classical products, won the national patents, and at the same time show the products on many International Exhibition of different countries.
- Shenchen Pump with good international reputation, signed Credit Guarantee Agreement with alibaba, and alibaba supply international trade security payment.
- Baoding peristaltic pump Engineering Technology Research Center officially established in shenchen pump, this is the only authority research institution in Baoding peristaltic pump industry.
- "Shenchen Brand" has become "China Peristaltic Pump Leading Brand", which can replace imported brand and conform to European standards, and has been **exported to 63 countries around the world.**

SHENCHEN

► Baoding Shenchen Precision Pump Co., Ltd

Baoding Peristaltic Pump Engineering Technology Research Center

Baoding Shenchen Precision Pump Co., Ltd is a high-tech enterprise, specialized in R&D, manufacturing and sales of peristaltic pump, syringe pump and other fluid precision transfer equipment. With rich technical force, excellent technological process, well corporate reputation, outstanding product quality, Shenchen constantly research and develop advance technology product, provide whole fluid solution program for customers. Shenchen is the professional manufacturer who can completely provide pump head, complete machine, OEM and ODM product.

"Shenchen team" is a high-qualified, young and innovative team, pose-ss research engineer, application engineer and professional sales and service team. In 2014, Baoding city government highly appreciated Shenchen's outstanding contribution in peristaltic pump industry. Government and Shenchen Precision Pump Co., Ltd, co-establish "Baoding Peristaltic Pump Engineering Technology Research Center", this center is joint approval by technology bureau, financial bureau, Development and Reform Commission, also the only authority research organization in .peristaltic pump industry.

"Shenchen product " is widely used in research laboratories, bio-pharmaceuticals, food and beverage, fine chemical engineering, environmental.etc, product get CE certification and exported to Europe, America, Australia, Southeast Asia, in total 63 countries.

"Shenchen Brand" has become "China Peristaltic Pump Leading Brand", which can replace imported brand and conform to European standards.



► Customer Feedback, from USA, UK, Italy, France, Korea.etc.

发送时间: 2015-06-04 09:59:47

收件人: Laura-ShenChen Pump

抄送:

主题: Pump arrived

Dear Miss Laura,

We have received the pump already. We also have tested the pump, we think the pump is very good. We will try to approach our customers and hopefully we can generate more inquiry and PO.

Thank you very much.

收件人: Amy Niu * <amy@good-pump.com>

抄送:

时间: 2015年12月11日 (周五) 01:55

大小: 124KB

Dear Amy,

Yes pump are working great, my customer is very pleased.

Kind regards,

日期: 2015年11月12日 星期四 18:51:13

主题: Re: Shipping documents

Dear Laura, I have received and checked the items. Thanks a lot.

The pumps are all ok. I have tested them.

Kindest regards,

Juan

收件人: Amy Niu * <amy@good-pump.com>

抄送:

时间: 2015年12月9日 (周三) 22:25

大小: 213KB

Hi Amy,

Pump has worked amazingly well. We don't need more pumps right now but

I will contact you when we need.

BR

Exported to 63 countries



Mainly exported to: American, Britain, France, Germany, Italy, Spain, Canada, Australia, Russia, India...



1. Cosmetic Production Line
2. Industrial Field
3. Filling Production Line
4. Laboratory Application
5. Analytical Instrument Supporting
6. Filling Production Line



1 Company	
2 Honor	
3 Application	
4 Products	
◆ OCM Products.....	2
◆ Compact Peristaltic Pump.....	4
◆ Flow Rates Peristaltic Pump	6
Lab2015.....	6
LabV Series/V Series.....	7
V6 Series (servo motor)	9
External control function.....	10
TC50.....	11
◆ Dispensing Peristaltic Pump.....	12
LabF Series/F Series.....	12
F6 Series (servo motor)	14
IF3 Low Pulsation (servo motor)	16
◆ Filling system.....	17
Split Type CF600/CF600 Plus.....	17
Split Type CF350/CF350 Plus (servo motor)	18
Integrated Type DF600/DF600 Plus.....	19
Integrated Type DF600IV/DF350 (servo motor)	20
◆ Standard Peristaltic Pump.....	26
LabN Series.....	26
BT-N Series.....	27
N6 Series (servo motor)	28
◆ Basic Peristaltic Pump.....	29
LabM Series.....	29
BT-M Series.....	30
◆ Explosion Proof Peristaltic Pump.....	32
EXP600.....	32
QD600.....	33
◆ Medical Dialysis Pump-UC15.....	34
◆ Compact Pump Head-MiniPump Series.....	35
◆ Multichannel Pump Head.....	37
DG Series.....	37
MC Series.....	39
MD Series.....	40
MT12	41
◆ Easy Load Pump Head-YZ Series.....	42
◆ Fast Load Pump Head-DZ Series.....	44
◆ Low Pulse Pump Head-DY Series.....	45
◆ Standard Pump Head-SN Series.....	46
5 OEM Products.....	47
6 Peristaltic Pump Tubing and Accessory.....	56

Revolutionary breakthrough in fluid control area

Shenchen provide **[OCM]** customized fluid solution

OCM (Original Customized Manufacture) customized fluid solution, with peristaltic pump core technology. Two working modes: Independent working mode and logic working mode.



- Same volume filling, different volume filling and transferring fluid
- Different units can co-work according with preset program.
- Independent working mode and logic working mode can work cooperatively.
- Different filling line, different dispensing request, choose freely, barrier-free conversion.



- Application1 Reactor liquid adding, dispensing and washing
- Application2 Dispensing different reagents in genetic engineering
- Application3 Pigment accurate allocation in dyeing and printing industry
- Application4 Several times filling in filling process to prevent foaming

- Peristaltic pump core technology
- One to one customized service
- Self-define working mode
- Automatic mixed filling follow commands
- Fast response, high reliability
- 7 inch- True color - Industrial LCD screen



Website: www.good-pump.com

Tel: 0086-312-6780680

Fax: 0086-312-6780636

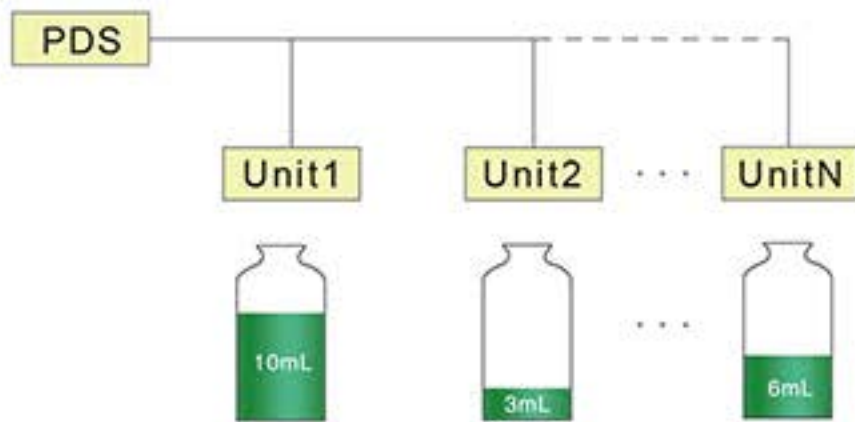


Diagram 1

In independent working mode, each unit working independently, can fill different volume.

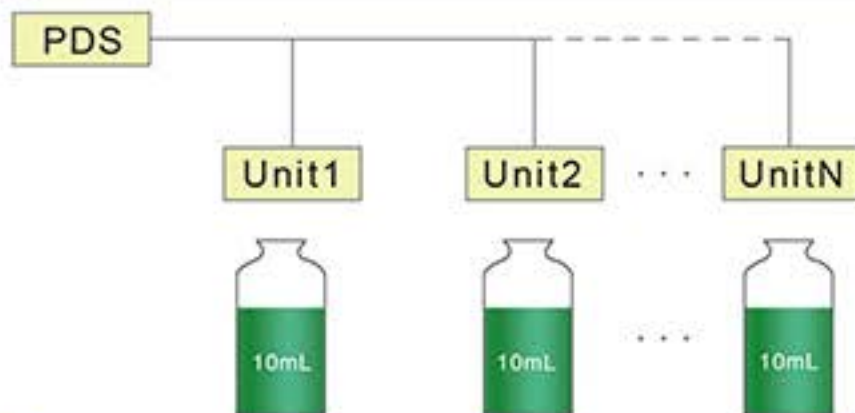


Diagram 2

In independent working mode, each unit can repeat filling same volume in same time, achieving normal dispensing function.

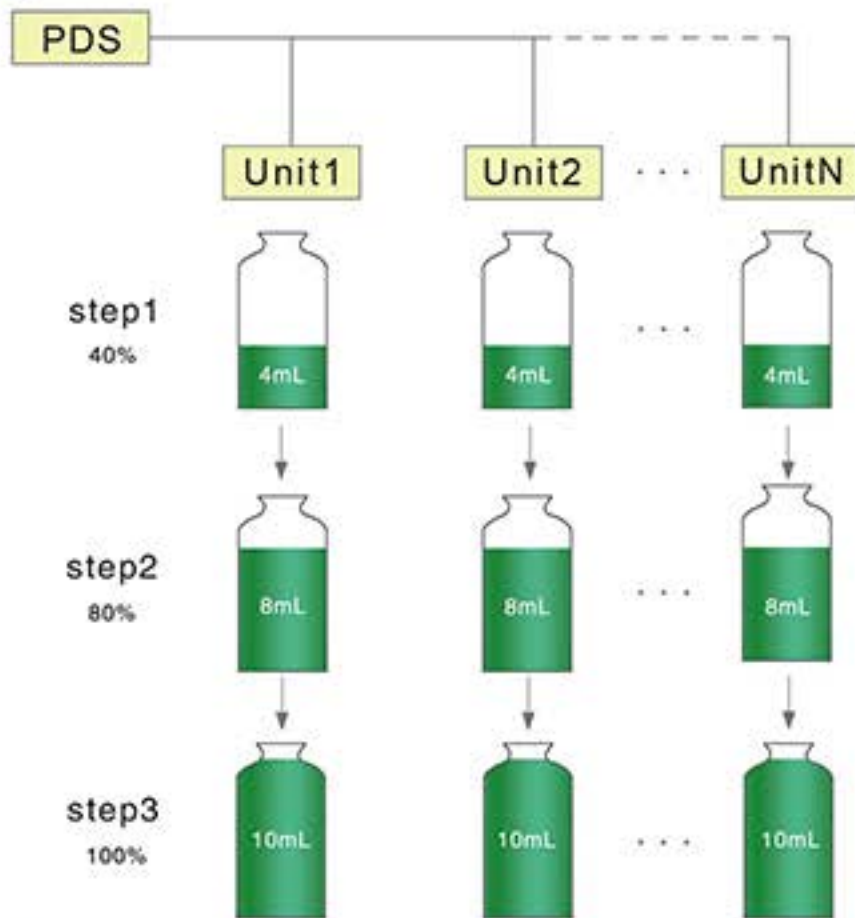


Diagram 3

In independent working mode, each unit can finish filling process with different filling volume. This function is suitable for prevent fluid splashing and foaming. This process can also repeat many times, achieving complicated dispensing function.

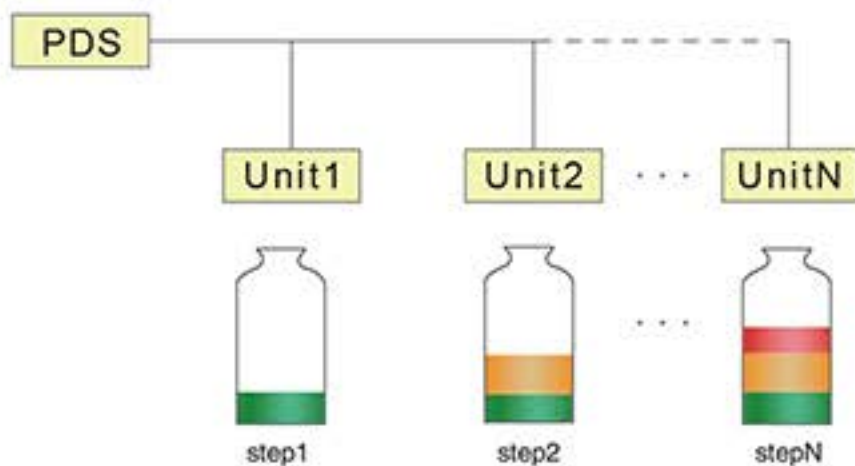


Diagram 4


In logic working mode, each unit working cooperatively. Filling or transferring with orders, also with different proportion adding, achieving automation mixing work.




Compact Peristaltic Pump SP-MiniPump



Installation Sizes Drawing (Unit: mm)

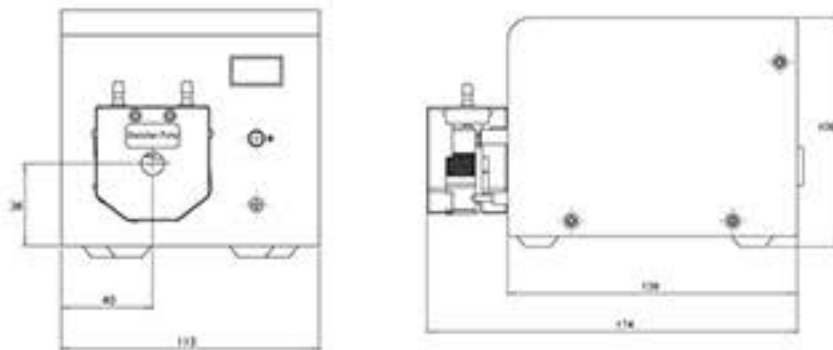

External Control Schematic Diagram

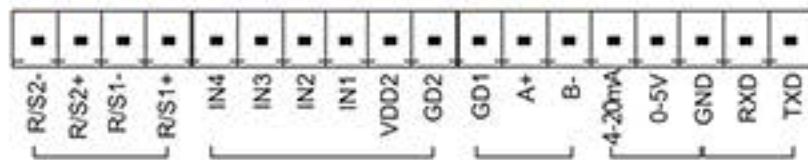

Features

- 304 stainless steel housing, resisting corrosion, space-saved.
- OLED screen, display the current motor speed and working status.
- Digital knob control speed, toggle switch control direction and start/stop.
- Various external control functions, support RS485 standard MODBUS protocol.

Model Number

- SP-MiniPump01
- SP-MiniPump02

Technical Specifications

Flow rate range	0.0024~190 mL/min
Speed range	0.1~300 rpm, reversible.
Speed resolution	0.1rpm
Speed control	Digital knob control speed
Motor type	Stepper motor
Display	OLED display
External control	Start/stop control (switch signal) Speed control 0-5V, 0-10V, 4-20mA for option(Analog signal)
Communication interface	RS485 support MODBUS protocol (RTU mode)
Power supply	24V/1A DC power
Drive dimension	136 × 113 × 108 mm (L × W × H)
Weight	1.64 kg
Power consumption	< 15 W
Condition temperature	0~40℃
Relative humidity	< 80%
IP rate	IP31

EM-MiniPump and SP-MiniPump speed/flow rates reference

Rotation speed (rpm)	Dual channel Flow rates (mL/min) and tubing (1 mm wall thickness)			
	ID 1 mm	ID 2 mm	ID 2.5 mm	ID 3 mm
5	0.25	0.90	1.28	1.78
10	0.49	1.84	2.59	3.52
15	0.74	2.76	3.85	5.30
25	1.25	4.57	6.43	8.80
40	2.00	7.31	10.25	14.02
50	2.52	9.14	12.78	17.47
75	3.76	13.64	19.04	26.12
100	5.03	18.17	25.42	35.34
200	10.12	36.45	51.27	71.44
250	12.63	45.75	64.34	88.35
300	15.01	54.63	76.84	108.39

Rotation speed (rpm)	Single channel Flow rates (mL/min) and tubing (1.6 mm wall thickness)			
	13# Tubing ID 0.8 mm	14# Tubing ID 1.6 mm	16# Tubing ID 3.1 mm	25# Tubing ID 4.8 mm
5	0.12	0.56	1.97	3.26
10	0.26	1.15	3.90	6.41
15	0.39	1.74	5.80	9.52
25	0.66	2.88	9.59	15.74
40	1.08	4.59	15.29	25.15
50	1.34	5.69	18.96	31.14
75	2.00	8.54	28.36	46.55
100	2.73	11.37	37.66	61.98
200	5.46	22.76	76.02	124.14
250	6.81	28.40	95.08	153.63
300	8.28	33.88	114.31	190.00



SHENCHEN Compact Peristaltic Pump

Compact Peristaltic Pump EM-MiniPump



Embedded design
-Single channel
EM-MiniPump01
(Install sample)


Features

- Embedded design, space-saved, much easier installation.
- OLED screen, display the current motor speed and working status.
- Membrane keypad control.


Model Number

- EM-MiniPump01
- EM-MiniPump02

Installation Sizes Drawing (Unit: mm)

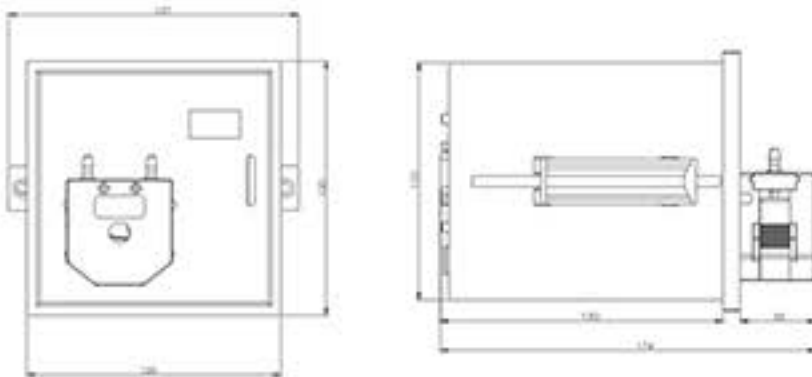

External Control Schematic Diagram

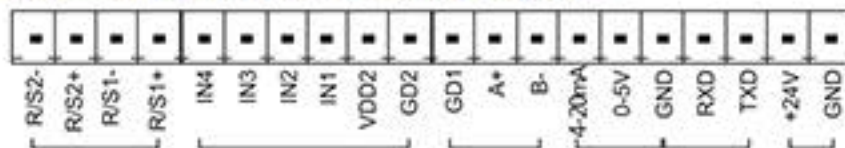

4 kinds of connect tubing way



Technical Specifications

Flow rate range	0.0024~190 mL/min
Speed range	0.1~300 rpm, reversible.
Speed resolution	0.1rpm
Control method	Membrane keypad
Motor type	Stepper motor
Display	OLED display
External control	Start/stop, direction control (switch signal) 0-5V, 0-10V, 4-20mA for option(Analog signal)
Communication interface	RS485 support MODBUS protocol(RTU mode)
Power supply	24V DC power
Drive dimension	110 × 120 × 120mm (L×W×H)
Weight	0.8 kg
Power consumption	< 15 W
Condition temperature	0~40℃
Relative humidity	< 80%
IP rate	IP31

Product Composition and Flow Rate Range

Drive	Motor type	Pump Head	Tubing size	Rotation speed	Flow rates range
EM-MiniPump	Stepper motor	MiniPump01 (4 rollers single channel)	1.6mm Wall thickness tubing: 13#, 14#, 16#, 25#	0.1~300 (rpm)	0.0024~190.00(mL/min)
			1mm Wall thickness tubing:Tubing ID 1mm, 2mm, 2.5mm, 3mm		
		MiniPump02 (4 rollers dual channel)	1mm Wall thickness tubing: Tubing ID 1mm, 2mm, 2.5mm, 3mm		0.005~108.39(mL/min)

Lab2015

(Glory Edition)



LAB
Intelligent
model



- Touch screen control
- Animation display
- ABS engineering plastic housing

Lab2015 Product Composition and Flow rate Range



Drive	Speed Range(rpm)	Pump Head	Tubing Size	Flow Rate(mL/min)
Lab2015	0.1-100rpm	YZ1515x	13#、14#、19#、16# 25#、17#、18#	0.007~380
		YZ2515x	15#、24#	0.17~290
		MC1-MC8(10)	Inner diameter: 0.13~3.17mm Wall thickness: 0.8~1mm	0.000166~32
		MC1-MC8(6)		0.000185~42



SHENCHEN Flow Rates Peristaltic Pump

Flow Rates Peristaltic Pump

LabV1, LabV3, LabV6



Touch screen control, Dynamic display

ABS engineering plastic housing

LAB Intelligent Type



Features


- Industrial grade 4.3" true color LCD screen, touch screen control, both Chinese and English interface available.
- Dynamic display transferring status. Flow rate data, setting parameters and system configuration display in the same screen.
- 3 Kinds of working mode: fixed volume metering, fixed time and volume, timer start and stop, meet different transferring and dispensing request.
- Intelligent calibration function and online micro adjusting function.

Model Number

- LabV1
- LabV3
- LabV6

Typical Application

- Laboratory
- Industrial equipment supporting

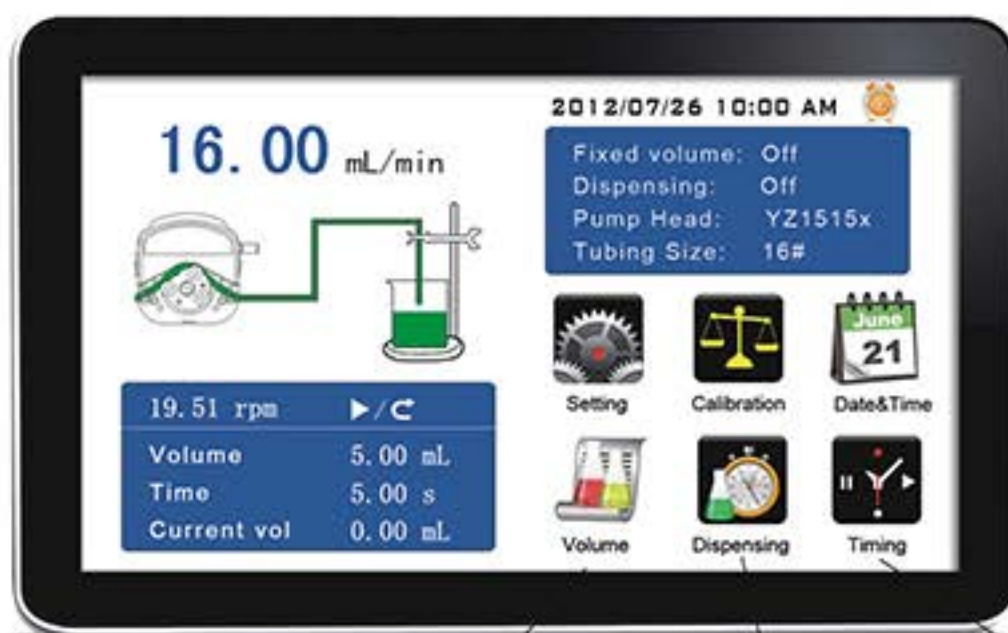

Flow Setting Interface



External Control Speed Setting Interface



Fixed Volume Measurement Setting Interface



Fixed Time and Volume Setting Interface



Timer Start and Stop Setting Interface



Flow Rates Peristaltic Pump V1,V3,V6



Touch screen
control, Dynamic
display

304
stainless steel
housing

LAB
Intelligent Type



Suitable Pump Head



Multichannel
Pump Head
(MC1-MC12)



Multichannel
Pump Head
(DG1-DG12)



Easy Load
Pump Head
(YZ1515x)



Easy Load
Pump Head
(YZ2515x)


Model Number

■ V1 ■ V3 ■ V6

Typical Application

■ Laboratory ■ Industrial equipment

External Control Schematic Diagram

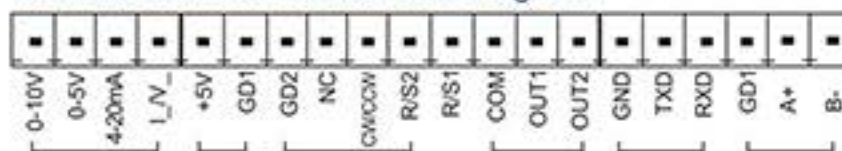

Technical Specifications

Flow rate range	LabV1/ V1:0.000166~570 mL/min LabV3/ V3:0.000166~1330 mL/min LabV6/ V6:0.000166~2280 mL/min	Start/stop, direction signal	Passive switch signal,such as foot pedal switch Active switch signal: 5V,12V and 24V for option
Speed range	0.1~600 rpm	Communication interface	RS232,RS485 support MODBUS protocol (RTU mode)
Speed resolution	0.01 rpm	Output interface	Output motor working status (Open-Collector output)
Flow rate resolution	0.01 μ l	Power supply	AC 220V $\pm$ 10% 50Hz/60Hz (standard) AC 110V $\pm$ 10% 50Hz/60Hz (optional)
Flow rate accuracy	< $\pm$ 0.5%	Drive dimension (L x W x H)	LabV Series: 261.4 x 157.3 x 236.9 mm V Series: 219.6 x 152 x 211.6 mm
Back suction angle	0~360°	Drive weight	LabV Series: 4.40 kg V Series: 4.20 kg
Outlet pressure	0.1Mpa (0.86~1.0mm wall thickness tubing) 0.1~0.27Mpa (1.6~2.4mm wall thickness tubing)	Power consumption	< 50W
Motor type	Stepper motor	Condition temperature	0~40°C
Display	Industrial grade 4.3" LCD color display	Relative humidity	< 80%
Control method	Touch screen and membrane keypad	IP rate	IP31
External speed control signal	0-5V,0-10V,4-20mA		

Product Composition and Flow Rate Range

Flow Rates Peristaltic Pump		Pump head & Flow rates (mL/min)			
		YZ1515x	YZ2515x	MC1-MC12(10)	MC1-MC12(6)
Drive & speed	Tubing	13', 14', 19', 16' 25', 17', 18'	15', 24'	Inner diameter: 0.13~3.17mm Wall thickness: 0.86~1mm	
LabV1/V1	0.1~150rpm	0.007~570	0.17~435	0.000166~49(suggest working rotating speed $\leq$ 150rpm)	0.000185~65(suggest working rotating speed $\leq$ 150rpm)
LabV3/V3	0.1~350rpm	0.007~1330	0.17~1015		
LabV6/V6	0.1~600rpm	0.007~2280	0.17~1740		



SHENCHEN Flow Rates Peristaltic Pump

Flow Rates Peristaltic Pump



Features

- Large flow rate, High precision, Intelligent control of liquid transferring.
- Servo motor drive, accurate control, strong driving force.
- 304 stainless steel housing, the first choice for high level industrial sites.

Model Number

- V6-3L/DZ25-3L
- V6-6L/DZ25-6L
- V6-12L/YZ35

Technical Specifications

Flow rate range	V6-3L: 0.211~3600 mL/min V6-6L: 0.3~6000 mL/min V6-12L: 0.00069~12 L/min	Start/stop,direction signal	Passive switch signal,such as foot pedal switch Active switch signal: 5V,12V and 24V for option
Speed range	0.1~600 rpm	Communication interface	RS232,RS485 support MODBUS protocol (RTU mode)
Speed resolution	0.01 rpm	Output interface	Output motor working status (Open-Collector output)
Flow rate resolution	0.01 mL	Power supply	AC 220V ± 10% 50Hz/60Hz (standard) AC 110V ± 10% 50Hz/60Hz (optional)
Flow rate accuracy	< ± 0.5%	Drive dimension (L × W × H)	V6-3L: 223 × 152 × 200mm ; V6-6L: 283 × 192 × 218mm V6-12L: 302 × 222 × 275mm
Back suction angle	0~360°	Drive weight	V6-3L: 5.02kg; V6-6L: 7.85kg; V6-12L: 13.14kg;
Outlet pressure	0.3 Mpa	Power consumption	V6-3L: < 80W ; V6-6L: < 180W ; V6-12L: < 300W
Motor type	Servo motor	Environment temperature	0~40℃
Display	Industrial grade 4.3" color LCD display	Relative humidity	< 80%
Control method	Touch screen and membrane keypad	IP rate	IP31
External speed control signal	0~5V,0~10V,4~20mA		

Product Composition and Flow Rate Range


Drive	Motor Type	Pump Head	Tubing Size	Speed Range(rpm)	Flow Rate Range(mL/min)
V6-3L	Servo Motor	DZ25-3L	15", 24", 35", 36"	0.1~600	0.211~3600
V6-6L		DZ25-6L	15", 24", 35", 36"		0.3~6000
V6-12L		YZ35	26", 73", 82"		0.69~12000



V Series peristaltic pump has various external control interface, can meet different equipment supporting requirements.

External control interface definition


1. Analog signal input port: Choose the external control speed signal from the external control setting interface. Turn on the external control speed function, control the motor speed range through the analog signal.
2. Internal isolation 5VDC output.
3. External control start/stop、direction signal input port: active signal input.
4. R/S1 external control start/stop signal input port: passive signal input.
5. Motor running status output port: output current running status of the motor.
6. RS232 communication interface: choose RS232 in communication setting interface, this port is effective.
7. RS485 communication interface: choose RS485 in communication setting interface, this port is effective.

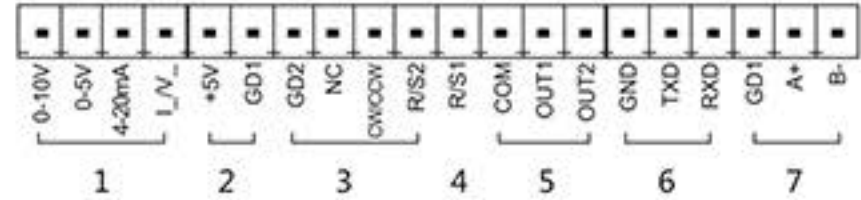

V Series peristaltic pump external control setting interface.

Choose external control start/stop、direction signal

1. Can respectively set external control start/stop or reversing switch whether effective or not.
2. Can choose control mode according to requirement: pulse mode or level mode.
3. Can set the foot switch whether effective or not.
4. Can choose switch value signal 's high level values according to customers' main control equipment: 5V(default)、12V、24V.

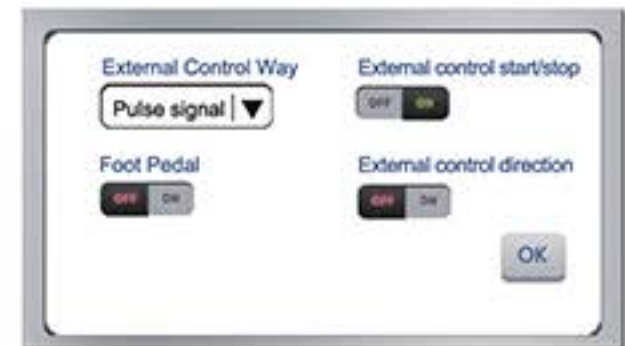


Programmable external control speed setting

1. Can choose analog signal source according to requirements: 0-5V, 0-10V, 4-20mA.
2. When the signal source is chosen, can choose maximum speed and minimum speed which corresponding to the signal source's maximum value and minimum value according to actual demand, to reach userdefined rotate speed range purpose.
3. Working speed limit is at the situation that the linear relation of analog quantity signal and rotating speed keeps invariant, set peristaltic pump's maximum working rotate speed. This setting can avoid production accident caused by sudden change of transfer fluid amount with the sudden change of external analog signal.




Communication Setting

1. Support RS232 and RS485 interface, can be chosen in software interface.
2. Support various Baud rate : 2400bps, 4800bps, 9600bps, 19200bps.
3. Can choose standard modbus communication protocol or Shencheng communication protocol to control, modbus communication protocol is suitable to match with industrial site HMI, PLC or other upper computer which support standard modbus communication protocol; Shencheng communication protocol is suitable to match with Singlechip micro computer or upper computer on computer programme.

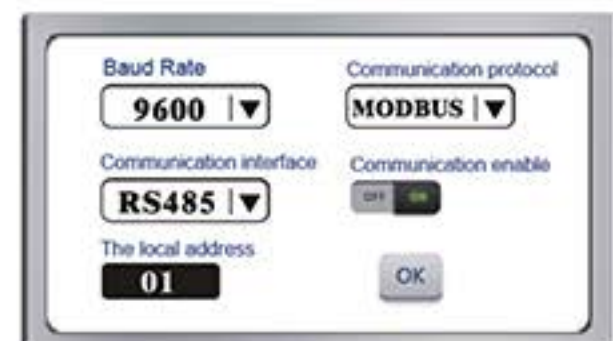

Installation Sizes Drawing (Unit: mm)

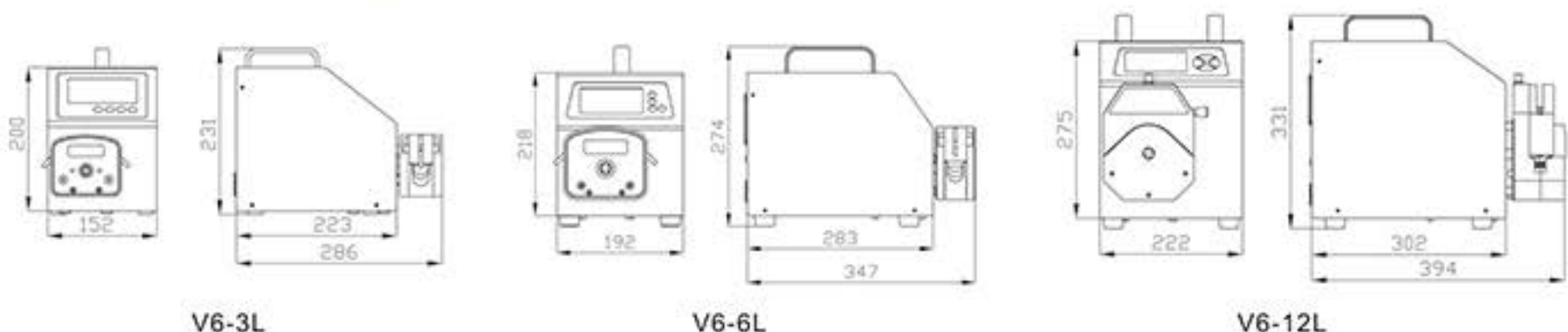



SHENCHEN

Flow Rates Peristaltic Pump

Flow Rates Peristaltic Pump TC50

304
stainless steel
housing



Special for flow injection analyzer

Features

- Special for flow injection analyzer.
- Overall block design to improve the accuracy between channels.
- The first choice for multichannel, high precision, and micro flow rate transferring.


Model Number ■ TC50



Technical Specifications

Flow rate range	0.0016~10.5 mL/min
Speed range	0.1~50 rpm
Speed resolution	0.01rpm
Flow rate resolution	0.01 μ l
Flow errors between channels	<3%
Back suction angle	0~360°
Outlet pressure	0.15 Mpa
Motor type	Servo motor
Display	Industrial grade 4.3" color LCD display
Control method	Touch screen and membrane keypad
External Speed control signal	0~5V, 0~10V, 0~20mA for option
Start/stop, direction signal	Passive switch signal, such as: foot pedal switch Active switch signal: 5V, 12V, 24V for option
Communication interface	RS232, RS485 support MODBUS protocol(RTU mode)
Output interface	Output motor working status (Open-Collector output)
Power supply	AC 220V $\pm$ 10% 50Hz/60Hz (Standard) AC 110V $\pm$ 10% 50Hz/60Hz (Optional)
Drive dimension	345 $\times$ 150 $\times$ 217mm (L $\times$ W $\times$ H)
Drive weight	8.02 kg
Power consumption	< 50 W
Condition temperature	0~40°C
Relative humidity	< 80%
IP rate	IP31

Installation Sizes Drawing (Unit: mm)

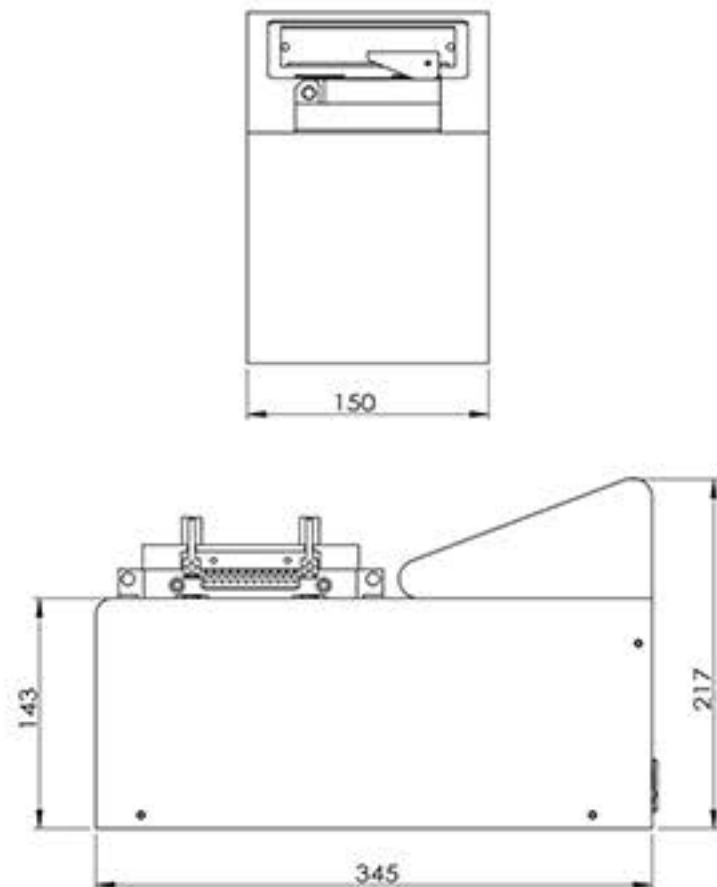

Product Composition and Flow Rate Range

Drive	Motor Type	Pump Head	Tubing Size	Speed Range(rpm)	Flow Rate Range(mL/min)
TC50	Servo Motor	MT12	0.51 $\times$ 0.9	0.1~50	0.0016~0.84
			1.02 $\times$ 0.86		0.005~2.5
			2.06 $\times$ 0.86		0.015~7.5
			2.54 $\times$ 0.86		0.020~10.5



Dispensing Peristaltic Pump

LabF1, LabF3, LabF6



Features

- Imported 4.3" industrial grade color display with touch screen control.
- Can preset dispensing volume, dispensing time, pause time and copy numbers.
- With intelligent calibration function and online micro adjusting function.
- The pump can store 60 commonly used filling modes.
- Back suction angle setting, avoid liquid drop off when the pump stops working.

Model Number

- LabF1
- LabF3
- LabF6



① Filling Countersunk ② Filling Nozzle ③ Foot pedal switch

Typical Application

- Medicine and chemical dispensing, such as oral liquid, diagnostic reagents.
- Cosmetic dispensing, such as perfume, lip gloss.



Dispensing volume interface



Dispensing time and pause time interface



Copy numbers interface



SHENCHEN Dispensing Peristaltic Pump

Dispensing Peristaltic Pump F1,F3,F6



Touch screen
control, Dynamic
display

304
stainless steel
housing



Model Number


■ F1 ■ F3 ■ F6

Technical Specifications

Flow rate range	LabF1/F1: 0.00166~570mL/min LabF3/F3: 0.00166~1330mL/min LabF6/F6: 0.00166~2280mL/min
Speed range	1~600 rpm
Speed resolution	0.01 rpm
Dispensing volume range	0.1~9999.99 mL
Dispensing volume resolution	0.01 mL
Dispensing time	0.5~9999.99 s
Pause time	0.5~9999.99 s
Time resolution	0.01 s
Copy numbers	1~9999 numbers, setting '0' means unlimited
Back suction angle	0~360°
Dispensing accuracy	< ± 0.5%
Motor type	Stepper motor
Display	Industrial grade 4.3" color LCD display
Control method	Touch screen and membrane keypad



External Control Schematic Diagram

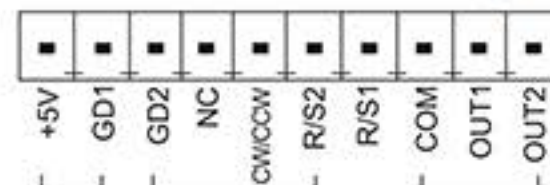

Product Composition and Flow Rate Range

Dispensing Peristaltic Pump		Pump head & Flow rates (mL/min)			
		YZ1515x	YZ2515x	MC1~MC12(10)	MC1~MC12(6)
Drive & speed	Tubing	13', 14', 19', 16', 25', 17', 18'	15', 24'	Inner diameter: 0.13~3.17mm Wall thickness: 0.86~1mm	
LabF1/F1	1~150rpm	0.07~570	1.7~435	0.00166~49(suggest working rotating speed ≤ 150rpm)	0.00185~65(suggest working rotating speed ≤ 150rpm)
LabF3/F3	1~350rpm	0.07~1330	1.7~1015		
LabF6/F6	1~600rpm	0.07~2280	1.7~1740		



Dispensing Peristaltic Pump



Features

- Suitable for large flow rate, high efficiency, high precision filling.
- Servo motor drive, accurate control, strong driving force.
- 304 stainless steel drive housing can support the filling line.

Model Number

- F6-3L/DZ25-3L
- F6-6L/DZ25-6L
- F6-12L/YZ35

Technical Specifications

Flow rate range	F6-3L: 2.11~3600 mL/min F6-6L: 3.0~6000 mL/min F6-12L: 0.0069~12 L/min	Display	Industrial grade 4.3" color LCD display
Speed range	1~600 rpm	Control method	Touch screen and membrane keypad
Speed resolution	0.01 rpm	Start/stop, direction signal	Passive switch signal, such as foot pedal switch Active switch signal: 5V, 12V and 24V for option
Dispensing volume range	0.1~9999.99 mL	Output interface	Output motor working status (Open-Collector output)
Dispensing volume resolution	0.01 mL	Power supply	AC 220V±10% 50Hz/60Hz (Standard) AC 110V±10% 50Hz/60Hz (Optional)
Dispensing time	0.5~9999.99 s	Drive dimension	F6-3L: 223×152×200mm ; F6-6L: 283×192×218mm F6-12L: 302×222×275mm
Pause time	0.5~9999.99 s	Drive weight	F6-3L: 5.02kg; F6-6L: 7.85kg; F6-12L: 13.14kg;
Time resolution	0.01 s	Power consumption	F6-3L: <80W ; F6-6L: <180W ; F6-12L: <300W
Copy numbers	1~9999 times, setting '0' means unlimited	Condition temperature	0~40°C
Back suction angle	0~360°	Relative humidity	< 80%
Outlet pressure	0.3 Mpa	IP rate	IP31
Dispensing accuracy	< ± 0.5%		
Motor type	Servo motor		

Product Composition and Flow Rate Range

Drive	Motor Type	Pump Head	Tubing Size	Speed Range(rpm)	Flow Rate Range(mL/min)
F6-3L	Servo Motor	DZ25-3L	15', 24', 35', 36'	1~600	2.11~3600
F6-6L		DZ25-6L	15', 24', 35', 36'		3.0~6000
F6-12L		YZ35	26', 73', 82'		6.9~12000



SHENCHEN

Dispensing Peristaltic Pump

LabF6/F6 Filling Volume Reference Parameter(Media is water)

Drive	Pump head	Tubing	Filling volume (mL)	Filling time(s)	Filling Accuracy (± %)	Output(pcs/min)	Motor Speed(rpm)
LabF6/F6	YZ1515x YZ2515x	13#	0.1	0.5	± 5ul	40	204.083
		13#	0.3	0.7	1.5	35	426.251
		13#	0.5	1	0.8	30	516.081
		13#	1	2	0.5	20	517.152
		14#	2	1	1	30	446.724
		14#	3	1.5	0.8	24	446.479
		19#	5	1.2	1	27	454.919
		16#	7	1	0.5	30	457.705
		25#/15#	10	1	1	30	303.426
		25#/15#	15	1	0.8	30	461.273
		25#/15#	20	1.2	0.5	27	518.945
		17#/24#	30	1.2	0.8	27	462.725
		17#/24#	50	2	0.5	20	461.595
		18#	80	2.5	0.5	17	427.274
		18#	100	3	0.5	15	446.583
F6-3L	DZ25-3L	15#	16	1	0.5	30	443.540
		24#	30	1.2	1.0	27	454.877
		35#	150	4	0.6	12	447.940
		36#	200	4	0.6	12	481.802
F6-6L	DZ25-6L	15#	80	4	0.4	12	396.800
		24#	150	4	0.4	12	440.700
		35#	200	3.2	0.5	14	439.540
		36#	300	3.5	0.5	13	473.208
F6-12L	YZ35	26#	150	3	0.5	15	423.254
		73#	300	3	0.5	15	457.805
		82#	500	3	0.5	15	458.451

F6 Series Filling Accuracy Reference Parameter(Media is water)

Drive	Pump head	Tubing	Filling volume	Filling time	Actual Filling Volume(mL)										Filling Accuracy	
F6-3L	DZ25-3L	15'	16mL	1s	15.97	15.97	16.00	16.06	15.97	15.92	15.96	16.02	16.01	15.97	-0.50%	0.50%
					15.98	16.00	16.04	16.08	15.92	15.95	16.01	16.01	15.97	15.97		
		24'	30mL	1.2s	29.98	30.12	30.20	30.29	30.06	29.94	30.11	30.18	30.02	29.99	-0.20%	1.00%
					30.23	30.24	30.12	29.99	30.12	30.23	30.20	30.30	30.19	30.11		
		35'	150mL	4s	149.50	149.50	150.10	150.00	149.80	149.00	149.70	149.50	150.10	149.60	-0.15%	0.60%
					149.00	149.90	149.70	149.10	149.50	150.00	150.10	150.10	149.60	149.40		
		36'	200mL	4s	199.70	200.30	200.10	200.30	200.00	200.00	200.70	200.80	200.70	200.40	-0.67%	0.07%
					200.30	200.00	200.30	200.80	200.90	200.70	200.80	200.40	200.50	200.20		
F6-6L	DZ25-6L	15'	80mL	4s	80.2	80.0	80.3	80.2	80.1	80.1	80.2	80.1	80.2	80.1	-0.125%	0.375%
					80.1	79.9	80.1	79.9	80.0	80.1	19.9	80.1	80.1	80.0		
		24'	150mL	4s	149.9	150.2	150.0	150.1	150.0	150.3	150.0	150.1	150.3	150.0	-0.067%	0.200%
					149.9	150.1	150.0	150.0	150.1	150.0	150.1	149.9	150.2	150.2		
		35'	200mL	3.5s	199.7	199.9	200.2	200.2	200.2	199.6	199.8	199.8	200.1	199.9	-0.400%	0.100%
					199.5	199.2	199.5	199.9	199.5	199.3	199.8	199.3	199.6	199.9		
		36'	300mL	3.5s	300.4	300.1	300.3	300.3	300.4	300.1	300.3	300.4	300.6	300.2	-0.167%	0.200%
					300.2	300.0	300.1	300.2	300.2	299.9	299.5	300.5	300.5	300.6		
F6-12L	YZ35	26'	150mL	3s	149.3	150.3	149.6	150.2	150.5	150.1	150.5	150.1	150.5	150.0	-0.47%	0.33%
					150.4	149.9	150.4	150.0	150.3	149.8	150.3	149.9	150.4	150.0		
		73'	300mL	3s	299.9	299.8	300.0	300.1	300.4	300.2	300.4	300.4	300.6	300.5	-0.07%	0.27%
					300.8	300.6	300.5	300.5	300.5	300.6	300.7	300.5	300.6	300.6		
		82'	500mL	3s	500.0	497.8	499.7	498.3	498.4	499.5	498.2	499.7	499.2	498.4	-0.44%	0.38%
					499.5	498.3	501.5	501.9	500.0	501.8	501.1	500.5	501.6	499.6		



Low Pulsation Dispensing Peristaltic Pump IF3



Servo
Motor Drive

304
stainless steel
housing


INDUSTRIAL
Intelligent Type



Features

- Servo motor drive, high precision, low pulsation dispensing peristaltic pump.
- The low pulsation pump head is special for high precision filling. Through the phase difference between the two sets of rollers, make the fluid peaks and valleys complementary, then reduce the pulsation of the fluid.
- Adaptive pressure tubing space, extend the tubing lifetime effectively.
- Achieving high precision dispensing of micro flow rate.

Installation Sizes Drawing (Unit: mm)

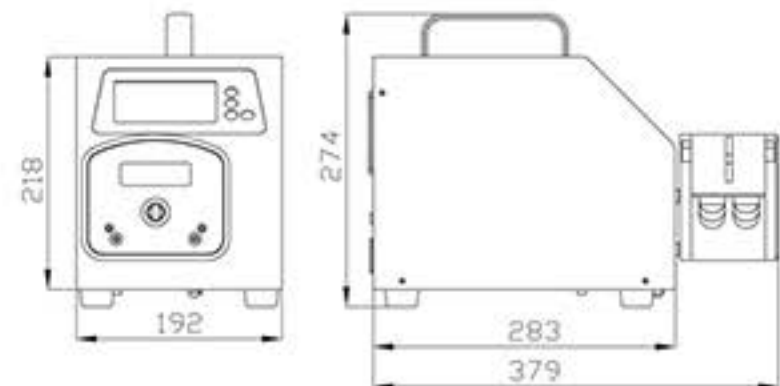

Model Number ■ IF3

Product Composition and Flow Rate Range

Drive	Motor Type	Pump head	Tubing	Speed range(rpm)	Flow rate range(mL/min)
IF3	Servo motor	DY15	13', 14', 19', 16', 25', 17', 18'	1-350	0.1-3337
		DY25	15', 24', 35', 36'		4.2-4340

Filling Volume Reference Parameter(Media is water)

Tubing	Filling volume	Filling time	Actual Filling Volume(mL)										Filling Accuracy	
			0.5003	0.5010	0.4982	0.5015	0.5015	0.5010	0.5014	0.5005	0.4981	0.4998	-0.37%	0.32%
14'	0.5mL	0.5s	0.5025	0.4983	0.5024	0.5021	0.5020	0.5004	0.5005	0.4997	0.4977	0.5008		
			1.0037	1.0000	1.0026	1.0003	1.0031	1.0015	1.0000	1.0007	1.0023	1.0035	-0.44%	0.42%
14'	1mL	1s	1.0002	1.0020	0.9968	1.0010	0.9986	1.0026	1.0030	1.0029	1.0023	1.0017		
			2.0020	2.0060	2.0039	2.0037	2.0047	1.9944	2.0062	2.0040	1.9951	2.0060	-0.58%	0.60%
14'	2mL	1s	2.0036	1.9940	2.0024	2.0086	1.9946	2.0005	1.9951	2.0068	1.9971	1.9965		
			3.0058	2.9935	3.0092	3.0008	2.9927	3.0072	3.0046	2.9914	3.0040	3.0000	-0.33%	0.31%
19'	3mL	1s	2.9902	3.0048	3.0008	2.9919		3.0012	3.0084	2.9915	2.9990	3.0051		
			5.0215	4.9982	5.0145	5.0038	4.9864	5.0183	4.9962	5.0227	5.0158	4.9806	-0.39%	0.50%
16'	5mL	1s	5.0248	5.0070	4.9861	5.0108	4.9995	5.0080	5.0044	4.9868	5.0231	4.9977		
			10.02	10.04	10.00	9.98	10.02	10.03	10.04	10.04	10.01	9.99	-0.20%	0.40%
25'	10mL	1s	10.00	10.02	10.02	10.04	9.99	9.98	10.00	10.04	10.02	10.03		
			15.02	14.97	14.97	14.94	15.00	15.00	14.94	14.95	14.96	14.99	-0.47%	0.13%
25'	15mL	1s	14.96	14.95	14.96	14.99	14.98	14.94	14.93	14.93	14.98	14.96		
17'	30mL	1s	29.92	30.01	29.99	29.91	29.83	29.86	29.91	29.91	29.91	29.89	-0.57%	0.20%
			29.88	29.96	30.03	30.06	30.02	30.02	30.06	29.96	29.96	29.83		
18'	100mL	2.5s	99.80	99.90	99.80	99.90	100.00	100.10	100.10	100.00	100.10	100.00	-0.20%	0.30%
			99.80	100.00	100.10	100.20	99.90	100.10	100.00	100.30	100.10	100.10		
15'	10mL	1s	9.97	10.01	9.96	9.98	10.01	9.96	9.97	9.99	9.96	9.96	-0.60%	0.10%
			10.00	9.94	9.97	10.00	9.96	9.97	10.00	10.00	10.01	9.97		
24'	30mL	1s	29.90	29.80	30.10	29.90	29.90	30.00	29.90	29.90	30.00	29.90	-0.67%	0.33%
			29.90	29.90	29.90	29.90	29.90	29.80	30.00	29.90	29.90	30.00		
35'	70mL	1.2s	69.80	69.90	69.90	69.90	69.70	70.10	70.00	69.70	69.90	69.90	-0.43%	0.43%
			69.90	70.10	70.00	69.90	70.30	69.70	70.10	69.70	69.70	69.90		
36'	100mL	2s	99.80	99.90	99.80	99.80	100.00	99.90	99.90	99.90	99.90	99.90	-0.20%	0.10%
			99.80	99.90	99.90	99.90	100.00	100.10	99.80	100.10	99.90	99.90		



Split Type Filling System

CF600

Touch screen control
dynamic display



Features

- Split design, adding or deleting filling units freely.
- One controller can control 1–16 filling units.
- Each filling units can receive stop filling signal when bottle absent.
- You can check the working status of the filling units from controller.


Model Number

- CF600
- CF600 Plus
- CF600 II
- CF600 II Plus


Typical Application

- Laboratory dispensing
- Industry filling



External Control Schematic Diagram

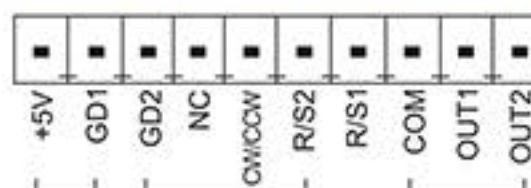


Installation Sizes Drawing (Unit: mm)

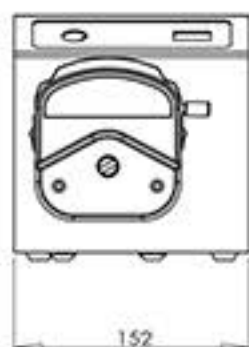

Drive Dimension Drawing

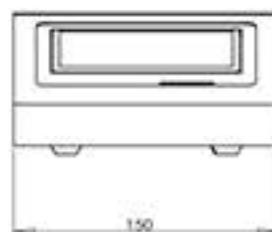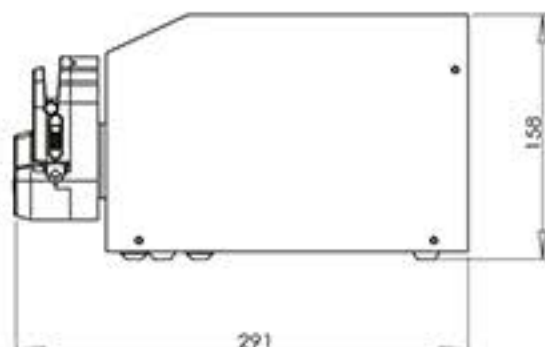

Controller Dimension Drawing

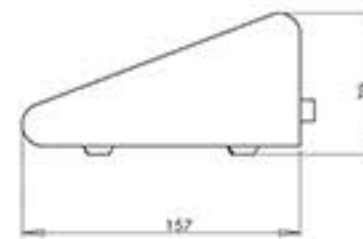



Split Type Filling system

CF350



Features

- Servo motor drive, precision control
- Split type design, space-saving, easy installation.
- Use low pulsation pump head to reduce the pulsation of the fluid effectively.
- Suitable for micro volume high precision filling.

Model Number



- CF350
- CF350 Plus
- CF600III
- CF600 PlusIII
- CF600IV
- CF600 PlusIV





Control Unit

Filling Unit

Installation Sizes Drawing (Unit: mm)

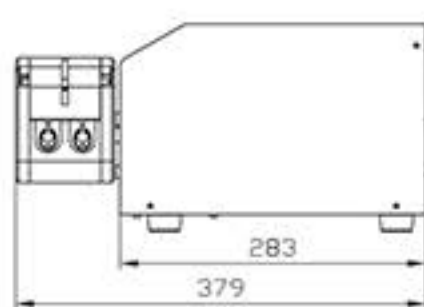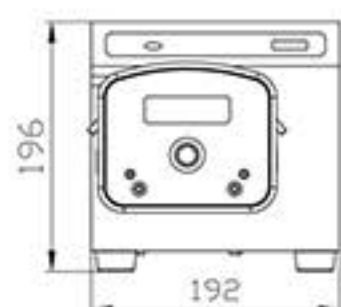

Drive Dimension Drawing

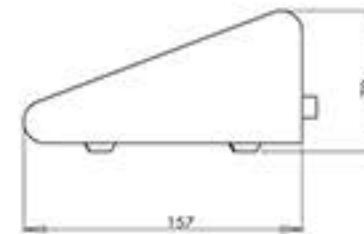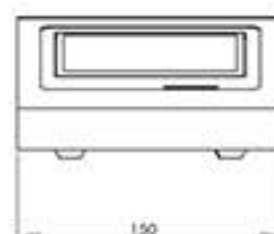

Controller Dimension Drawing

The pump type suffixed with "Plus" support MODBUS communication protocol.

Drive	Pump head	Flow rate range (mL/min)	Motor type	Drive dimension drawing (L x W x H)	Power consumption	Drive weight
CF600 CF600Plus	YZ1515x	0.07~2280	Stepper motor	237.4 x 152 x 158mm	Each unit<50W	4.2kg
	YZ2515x	1.7~1740				
CF600 II CF600Plus II	DZ25-3L	2.11~3600	Stepper motor	237.4 x 152 x 158mm	Each unit<80W	4.95kg
CF600 III CF600Plus III	DZ25-6L	3.0~6000	Servo motor	283 x 192 x 196mm	Each unit<180W	7.8kg
CF350 CF350Plus	DY15	0.1~3337	Servo motor	283 x 192 x 196mm	Each unit<180W	7.8kg
	DY25	4.2~4340				
CF600 IV CF600PlusIV	YZ35	6.9~12000	Servo motor	310 x 228 x 248mm	Each unit<300W	11.9kg
	YZ35-PPS					



Integrated Filling System

DF600



Features

- Integrated design, one controller can control 4 groups (total 16) filling units.
- It can support the production filling line, store 60 commonly used filling modes.
- Each filling units can independently receive stop filling signal when bottle absent.

Model Number

- DF600
- DF600 II
- DF600 Plus
- DF600 Plus II

Typical Application

- Pharmaceutical, health product filling, diagnostic reagent dispensing.
- Food, beverage filling.
- Cosmetic filling.

DF600 Plus



Features

Connected with HMI, PLC etc, industrial control system, one controller can control 4 groups (total 16) filling units. It can display the details of 16 filling units in the HMI, widely used in industrial control application.



Dispensing volume interface



Dispensing time and pause time interface



Copy numbers interface



DF600IV(servo motor control)

Servo
Motor Drive

High
precision filling

Software
Copyright
shenchen



DF350(servo motor control)

Servo
Motor Drive

High
precision filling

Software
Copyright
shenchen



Product Introduction

DF servo control series is split type intelligent filling system with high precision and low pulsation driven by servo motor. It is composed by control unit and integrated type filling units, each group have 4 filling units, total 16 filling units. This system use imported 4.3" industrial grade true color display with touch screen control, can preset filling volume, filling time, pause time, copy numbers and back-suction angle, animation display working status, filling data, setting parameter, system configuration display at the same screen; with intelligent calibration and online micro adjust function. Can connect foot pedal and receive switch signal, realize long-distance control. With motor working status output signal, can let other equipment realtime monitor the current filling status of peristaltic pump. This filling system use servo motor, and high precision, low pulsation pump head to make the filling precision higher, up to 0.1%–0.6%, make the micro volume and big volume high precision filling come true.

Model Number


- DF350
- DF350 Plus
- DF600III
- DF600 PlusIII
- DF600IV
- DF600 PlusIV

Typical Application


- Pharmaceutical, health product filling, diagnostic reagent dispensing.
- Food, beverage filling.
- Cosmetic filling.



Installation Sizes Drawing (Unit: mm)

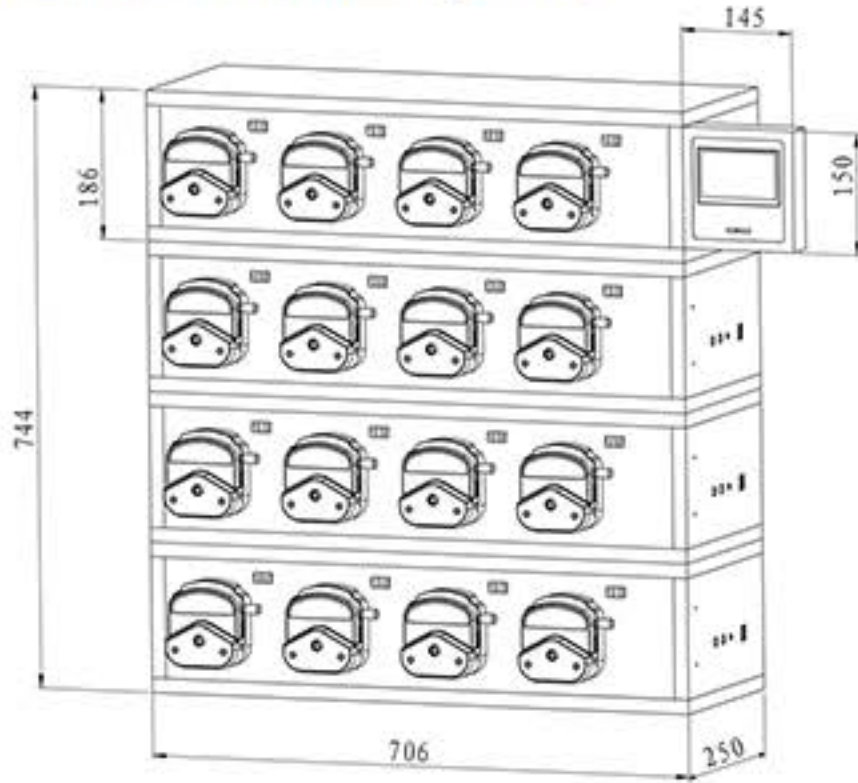


DF600 Dimension

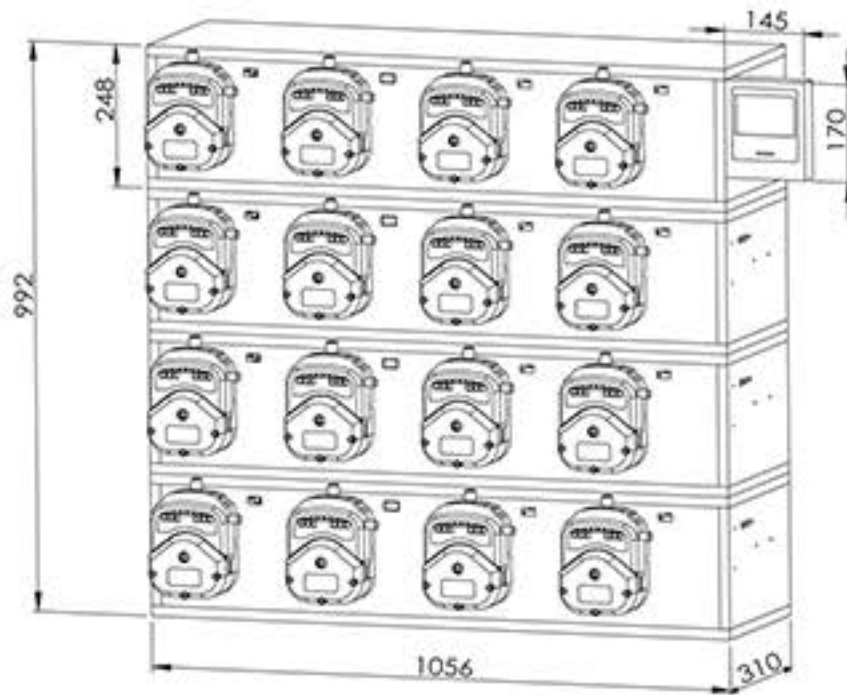


DF600IV Dimension

External Control Schematic Diagram


- DF600 has filling status output signal, can be monitored by other equipments.
- DF600 plus adopt standard MODBUS communication protocol, can be well connected with human-computer interface to achieve logic control.

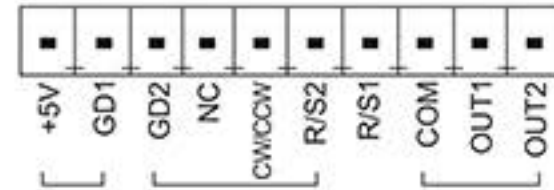

DF600 External Control Schematic Diagram

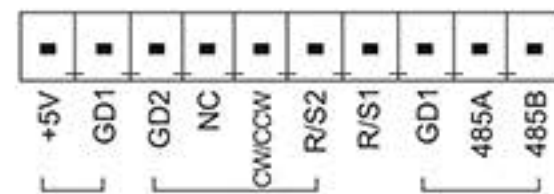

DF600plus External Control Schematic Diagram

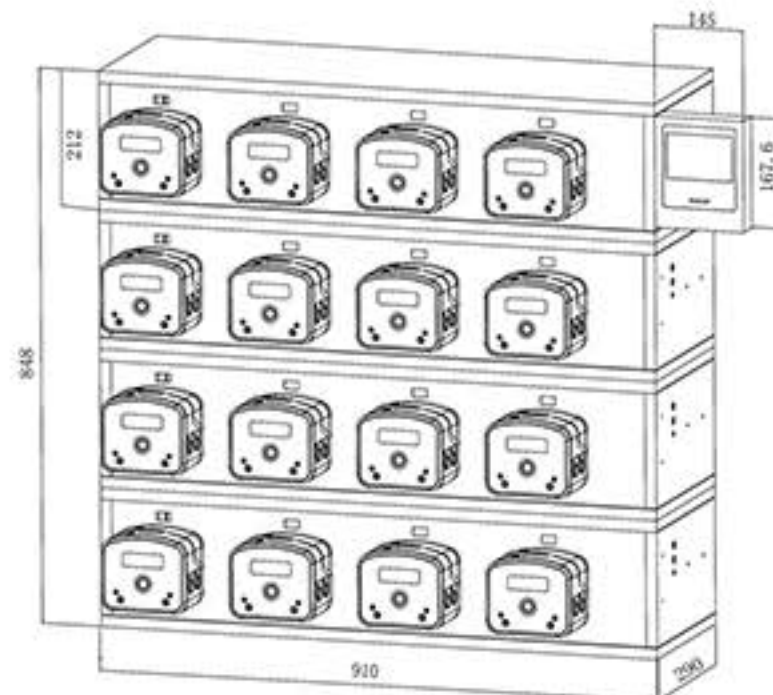

DF350 Dimension

The pump type suffixed with "Plus" support MODBUS communication protocol

Drive	Pump Head	Flow Rate Range (mL/min)	Motor Type	Drive Dimension (L x W x H)	Power Consumption	Drive Weight
DF600 DF600Plus	YZ1515x YZ2515x	0.07~2280 1.7~1740	Stepper motor	706 x 250 x 186mm	Each Unit<50W	21.8kg
DF600 II DF600Plus II	DZ25-3L	2.11~3600	Stepper motor	706 x 250 x 186mm	Each Unit<80W	21.8kg
DF600 III DF600Plus III	DZ25-6L	3.0~6000	Servo motor	910 x 290 x 212mm	Each Unit<180W	41.45kg
DF350 DF350Plus	DY15 DY25	0.1~3337 4.2~4340	Servo motor	910 x 290 x 212mm	Each Unit<180W	41.45kg
DF600 IV DF600Plus IV	YZ35 YZ35-PPS	6.9~12000	Servo motor	1056 x 310 x 248m	Each Unit<300W	58.6kg



Technical Specifications

Speed range	CF350/DF350	1-350 rpm	Filling accuracy	< ± 0.5%
	CF350 Plus/DF350 Plus	1-350 rpm	Display	4.3" industrial grade true color LCD screen
	CF600 III/DF600 III/DF600 IV	1-600 rpm	Control method	Touch screen and membrane keypad
	CF600 Plus III/DF600 Plus III/DF600 Plus IV	1-600 rpm	Start/stop, direction signal	Passive switch signal, such as foot pedal switch Active switch signal: 5V, 12V, 24V for option
Filling volume range	0.1-9999.99 mL		Output interface	Output motor working status (Open-Collector output)
Filling time range	0.5-9999.99 s		Power supply	AC 220V ± 10% 50Hz/60Hz (standard) AC 110V ± 10% 50Hz/60Hz (optional)
Pause time range	0.5-9999.99 s		Condition temperature	0-40°C
Filling volume resolution	0.01 mL		Relative humidity	< 80%
Time resolution	0.01 s		IP rate	IP31
Copy numbers	1-9999, '0' means unlimited			
Back suction angle	0-360°			

Product composition and flow rate range

Drive	Pump Head	Housing material	Tube clamp/ connector material	Tubing		(Motor speed 1-600rpm) Flow rate range mL/min
				Tubing size	ID x wall thickness	
DF600 CF600 DF600 Plus CF600 Plus	YZ1515x	PSF PPS	—	13'	0.8 × 1.6(mm)	0.07~42
				14'	1.6 × 1.6(mm)	0.27~162
				19'	2.4 × 1.6(mm)	0.55~330
				16'	3.1 × 1.6(mm)	0.82~492
				25'	4.8 × 1.6(mm)	1.7~1020
				17'	6.4 × 1.6(mm)	2.9~1740
	YZ2515x	PSF PPS	—	18'	7.9 × 1.6(mm)	3.8~2280
				15'	4.8 × 2.4(mm)	1.7~1020
DF600 II CF600 II DF600 Plus II CF600 Plus II	DZ25-3L	Aluminum alloy/PPS	PP	24'	6.4 × 2.4(mm)	2.9~1740
				15'	4.8 × 2.4(mm)	2.11~1264
				24'	6.4 × 2.4(mm)	3.85~2310
				35'	7.9 × 2.4(mm)	5.08~3050
DF600 III CF600 III DF600 Plus III CF600 Plus III	DZ25-6L	Aluminum alloy	PP	36'	9.6 × 2.4(mm)	6~3600
				15'	4.8 × 2.4(mm)	3~1800
				24'	6.4 × 2.4(mm)	5.5~3300
				35'	7.9 × 2.4(mm)	8~4800
DF600 IV CF600 IV DF600 Plus IV CF600 Plus IV	YZ35-PPS	PPS	PP	36'	9.6 × 2.4(mm)	10~6000
				26'	6.4 × 3.3(mm)	6.9~4200
				73'	9.6 × 3.3(mm)	12.3~7400
DF350 CF350 DF350 Plus CF350 Plus	DY15	Aluminum alloy	PP	82'	12.7 × 3.3(mm)	20~12000
				Tubing	ID x wall thickness	(Motor speed 1-350rpm) Flow rate range mL/min
	DY25	Aluminum alloy	PP	13'	0.8 × 1.6(mm)	0.1~48
				14'	1.6 × 1.6(mm)	0.6~223
				19'	2.4 × 1.6(mm)	1.3~448
				16'	3.1 × 1.6(mm)	2~723
				25'	4.8 × 1.6(mm)	4.7~1626
				17'	6.4 × 1.6(mm)	6.4~2230
				18'	7.9 × 1.6(mm)	9.5~3337
				15'	4.8 × 2.4(mm)	4.2~1480
				24'	6.4 × 2.4(mm)	7.6~2670
				35'	7.9 × 2.4(mm)	10~3600
				36'	9.6 × 2.4(mm)	12.4~4340



Filling Volume Reference (Media is water)

Drive	Pump Head	Tubing	Filling Volume (mL)	Filling Time (s)	Accuracy (± %)	Output (pcs/min)	Motor speed (rpm)
DF600 CF600 DF600 Plus CF600 Plus	YZ1515x YZ2515x	13#	0.1	0.5	± 5ul	40	204.083
		13#	0.3	0.7	1.5	35	426.251
		13#	0.5	1	0.8	30	516.081
		13#	1	2	0.5	20	517.152
		14#	2	1	1	30	446.724
		14#	3	1.5	0.8	24	446.479
		19#	5	1.2	1	27	454.919
		16#	7	1	0.5	30	457.705
		25#/15#	10	1	1	30	303.426
		25#/15#	15	1	0.8	30	461.273
		25#/15#	20	1.2	0.5	27	518.945
		17#/24#	30	1.2	0.8	27	462.725
		17#/24#	50	2	0.5	20	461.595
		18#	80	2.5	0.5	17	427.274
		18#	100	3	0.5	15	446.583
DF600 II CF600 II DF600 Plus II CF600 Plus II	DZ25-3L	15#	16	1	0.5	30	454.976
		24#	30	1.2	0.5	27	389.610
		35#	150	4	0.5	12	442.913
		36#	200	4	0.5	12	500.000
DF350 CF350 DF350 Plus CF350 Plus	DY15	14#	0.5	0.5	0.3	60	94.1915
		14#	1	1	0.4	30	94.1915
		14#	2	1	0.6	30	188.383
		19#	3	1	0.3	30	140.625
		16#	5	1	0.5	30	145.278
		25#	10	1	0.3	30	129.143
		25#	15	1	0.4	30	193.715
		17#	30	1	0.4	30	282.530
	DY25	18#	100	2.5	0.3	17	251.730
		15#	10	1	0.3	30	141.911
		26#	30	1	0.5	30	235.972
		35#	70	1.2	0.4	27	340.466
		36#	100	2	0.2	20	241.935
DF600 III CF600 III DF600 Plus III CF600 Plus III	DZ25-6L	15#	80	4	0.4	12	400.000
		24#	150	4	0.4	12	409.090
		35#	200	3.2	0.5	14	468.750
		36#	300	3.5	0.5	13	514.285
DF600 IV CF600 IV DF600 Plus IV CF600 Plus IV	YZ35-PPS	24#	150	3	0.5	15	428.570
		35#	300	3	0.5	15	486.499
		36#	500	3	0.5	15	500.000





Filling Volume Reference (Media is water)

	Tubing	Filling Volume	Filling Time	Actual Filling Volume(mL)										Filling Accuracy	
YZ1515x YZ2515x	13#	0.1mL	0.5s	0.1000	0.1017	0.1004	0.1000	0.1006	0.0986	0.1015	0.0988	0.1006	0.0991	-2.500%	1.700%
				0.0996	0.0975	0.0976	0.1006	0.0995	0.1015	0.0988	0.0991	0.0996	0.0993		
	13#	0.3mL	0.7s	0.3000	0.2979	0.3013	0.2985	0.3020	0.3030	0.3023	0.2993	0.2993	0.2974	-0.867%	1.033%
				0.3010	0.3030	0.2992	0.2979	0.3004	0.3019	0.3012	0.2978	0.3031	0.3022		
	13#	0.5mL	1s	0.5008	0.5034	0.4998	0.5020	0.5007	0.5038	0.5021	0.5023	0.5010	0.4995	-0.100%	0.760%
				0.5003	0.5008	0.5003	0.5011	0.5002	0.5037	0.5025	0.5005	0.5010	0.5000		
	13#	1mL	2s	1.0004	1.0000	0.9990	1.0011	1.0035	1.0014	1.0014	1.0014	0.9999	1.0015	-0.370%	0.490%
				1.0013	1.0013	1.0029	1.0049	1.0016	1.0019	1.0048	0.9963	1.0038	1.0040		
	14#	2mL	1s	2.0030	2.0030	2.0125	1.9861	2.0116	1.9987	1.9984	2.0143	1.9973	2.0110	-0.695%	0.765%
				2.0033	2.0153	1.9875	2.0100	2.0092	2.0024	2.0103	1.9915	2.0045	2.0133		
	14#	3mL	1.5s	3.0147	2.9861	3.0125	2.9897	3.0037	2.9905	3.0063	2.9850	3.0038	2.9917	-0.500%	0.493%
				3.0088	3.0000	3.0019	3.0138	2.9896	3.0148	2.9921	3.0128	2.9905	3.0059		
	19#	5mL	1.2s	4.9933	5.0069	4.9914	4.9540	5.0068	4.9957	4.9950	5.0029	4.9963	4.9571	-0.920%	0.500%
				4.9838	5.0208	5.0031	4.9687	5.0237	4.9841	5.0239	5.0032	4.9584	5.0250		
	16#	7mL	1s	7.0184	7.0197	7.0178	7.0120	7.0242	7.0132	7.0149	7.0213	7.0093	7.0215	-0.293%	0.346%
				6.9795	6.9834	6.9946	6.9975	6.9868	7.0133	7.0195	7.0121	7.0190	7.0125		
DZ25-3L	25# 15#	10mL	1s	9.9400	10.0300	10.0700	10.0200	10.0000	9.9800	9.9200	10.0400	10.0600	10.0400	-0.800%	0.700%
				10.0400	10.0500	10.0200	9.9700	9.9500	10.0200	10.0600	10.0400	10.0500	9.9600		
	25# 15#	15mL	1s	15.0000	15.0000	15.0300	14.9900	15.0400	15.0400	15.0300	15.0700	14.9000	15.0100	-0.667%	0.467%
				14.9500	15.0300	15.0300	15.0600	15.0300	14.9900	15.0300	14.9400	15.0100	15.0000		
	25# 15#	20mL	1.2s	19.9600	19.9400	20.0600	19.9900	19.9600	19.9200	19.9700	20.0400	20.0000	20.0000	-0.450%	0.500%
				19.9300	20.1000	20.0000	19.9100	20.0200	20.1000	20.0100	19.9700	19.9900	20.0900		
	17# 24#	30mL	1.2s	29.9100	29.9900	30.1800	30.1000	29.9700	29.9700	30.1800	29.8100	30.1800	29.9500	-0.633%	0.600%
				29.9400	29.9200	29.9300	30.0100	30.1800	30.0500	29.9600	29.8800	30.0100	30.1800		
	17# 24#	50mL	2s	50.0600	50.0400	49.9900	49.8500	50.0800	50.1100	49.9800	49.8700	49.9200	50.1200	-0.440%	0.240%
				49.9200	50.0000	49.9900	49.7800	49.8200	49.9100	50.0500	49.9300	49.8400	49.8400		
	18#	80mL	2.5s	79.9700	80.2200	80.2800	80.1100	80.0000	79.9000	79.9200	80.2900	80.1300	80.0300	-0.325%	0.363%
				79.7400	79.7800	79.8000	80.0100	80.1700	79.8500	79.7600	79.8500	80.0300	80.2200		
	18#	100mL	3s	99.8300	100.3500	99.9700	100.3700	100.0900	100.2600	100.0800	100.1500	100.1600	99.8900	-0.430%	0.370%
				100.0400	99.6600	100.0800	99.5700	100.0100	99.9900	99.7900	100.1400	99.7200	100.1700		
	35#	150mL	4s	149.4100	150.4300	149.7600	150.0400	150.4200	149.5500	150.3200	149.7500	149.6500	150.2900	-0.480%	0.287%
				149.4200	149.9400	150.2700	149.3400	149.8600	150.2500	149.2800	150.0300	149.9200	149.4700		
	36#	200mL	4s	199.4800	199.9000	199.6600	199.7200	199.8600	200.0200	199.2000	199.8800	199.8600	199.4800	-0.400%	0.300%
				199.9000	200.5400	200.6000	199.9600	199.9400	200.0200	199.9800	200.0600	199.9200	200.1800		





Filling Volume Reference (Media is water)

	Tubing	Filling Volume	Filling Time	Actual Filling Volume(mL)										Filling Accuracy	
DY15	14'	0.5mL	0.5s	0.5003	0.5010	0.4982	0.5015	0.5015	0.5010	0.5014	0.5005	0.4981	0.4998	-0.37%	0.32%
				0.5025	0.4983	0.5024	0.5021	0.5020	0.5004	0.5005	0.4997	0.4977	0.5008		
	14'	1mL	1s	1.0037	1.0000	1.0026	1.0003	1.0031	1.0015	1.0000	1.0007	1.0023	1.0035	-0.44%	0.42%
				1.0002	1.0020	0.9968	1.0010	0.9986	1.0026	1.0030	1.0029	1.0023	1.0017		
	14'	2mL	1s	2.0020	2.0060	2.0039	2.0037	2.0047	1.9944	2.0062	2.0040	1.9951	2.0060	-0.58%	0.60%
				2.0036	1.9940	2.0024	2.0086	1.9946	2.0005	1.9951	2.0068	1.9971	1.9965		
	19'	3mL	1s	3.0058	2.9935	3.0092	3.0008	2.9927	3.0072	3.0046	2.9914	3.0040	3.0000	-0.33%	0.31%
				2.9902	3.0048	3.0008	2.9919	2.9984	3.0012	3.0084	2.9915	2.9990	3.0051		
	16'	5mL	1s	5.0215	4.9982	5.0145	5.0038	4.9864	5.0183	4.9962	5.0227	5.0158	4.9806	-0.39%	0.50%
				5.0248	5.0070	4.9861	5.0108	4.9995	5.0080	5.0044	4.9868	5.0231	4.9977		
DY25	25'	10mL	1s	10.02	10.04	10.00	9.98	10.02	10.03	10.04	10.04	10.01	9.99	-0.20%	0.40%
				10.00	10.02	10.02	10.04	9.99	9.98	10.00	10.04	10.02	10.03		
	25'	15mL	1s	15.02	14.97	14.97	14.94	15.00	15.00	14.94	14.95	14.96	14.99	-0.47%	0.13%
				14.96	14.95	14.96	14.99	14.98	14.94	14.93	14.93	14.98	14.96		
	17'	30mL	1s	29.92	30.01	29.99	29.91	29.83	29.86	29.91	29.91	29.91	29.89	-0.57%	0.20%
				29.88	29.96	30.03	30.06	30.02	30.02	30.06	29.96	29.96	29.83		
	18'	100mL	2.5s	99.80	99.90	99.80	99.90	100.00	100.10	100.10	100.00	100.10	100.00	-0.20%	0.30%
				99.80	100.00	100.10	100.20	99.90	100.10	100.00	100.30	100.10	100.10		
	15'	10mL	1s	9.97	10.01	9.96	9.98	10.01	9.96	9.97	9.99	9.96	9.96	-0.60%	0.10%
				10.00	9.94	9.97	10.00	9.96	9.97	10.00	10.00	10.01	9.97		
DZ25-6L	24'	30mL	1s	29.90	29.80	30.10	29.90	29.90	30.00	29.90	29.90	30.00	29.90	-0.67%	0.33%
				29.90	29.90	29.90	29.90	29.90	29.80	30.00	29.90	29.90	30.00		
	35'	70mL	1.2s	69.80	69.90	69.90	69.90	69.70	70.10	70.00	69.70	69.90	69.90	-0.43%	0.43%
				69.90	70.10	70.00	69.90	70.30	69.70	70.10	69.70	69.70	69.90		
	36'	100mL	2s	99.80	99.90	99.80	99.80	100.00	99.90	99.90	99.90	99.90	99.90	-0.20%	0.10%
				99.80	99.90	99.90	99.90	100.00	100.10	99.80	100.10	99.90	99.90		
	15'	80mL	4s	80.2	80.0	80.3	80.2	80.1	80.1	80.2	80.1	80.2	80.1	-0.13%	0.38%
				80.1	79.9	80.1	79.9	80.0	80.1	19.9	80.1	80.1	80.0		
	24'	150mL	4s	149.9	150.2	150.0	150.1	150.0	150.3	150.0	150.1	150.3	150.0	-0.07%	0.20%
				149.9	150.1	150.0	150.0	150.1	150.0	150.1	149.9	150.2	150.2		
YZ35-PPS	35'	200mL	3.5s	199.7	199.9	200.2	200.2	200.2	199.6	199.8	199.8	200.1	199.9	-0.40%	0.10%
				199.5	199.2	199.5	199.9	199.5	199.3	199.8	199.3	199.6	199.9		
	36'	300mL	3.5s	300.4	300.1	300.3	300.3	300.4	300.1	300.3	300.4	300.6	300.2	-0.17%	0.20%
				300.2	300.0	300.1	300.2	300.2	299.9	299.5	300.5	300.5	300.6		
	26'	150mL	3s	149.3	150.3	149.6	150.2	150.5	150.1	150.5	150.1	150.5	150.0	-0.47%	0.33%
				150.4	149.9	150.4	150.0	150.3	149.8	150.3	149.9	150.4	150.0		
	73'	300mL	3s	299.9	299.8	300.0	300.1	300.4	300.2	300.4	300.4	300.6	300.5	-0.07%	0.27%
				300.8	300.6	300.5	300.5	300.5	300.6	300.7	300.5	300.6	300.6		
	82'	500mL	3s	500.0	497.8	499.7	498.3	498.4	499.5	498.2	499.7	499.2	498.4	-0.44%	0.38%
				499.5	498.3	501.5	501.9	500.0	501.8	501.1	500.5	501.6	499.6		





Standard Peristaltic Pump

LabN1, LabN3, LabN6



Features

- 3.2" blue LCD screen display, touch keypad control.
- Super silent drive, low vibration, low noise.
- Timing function, time range 0.1s-9999hours, can be used for simple dispensing function.

Model Number

- LabN1
- LabN3
- LabN6

Typical Application

- Special for university laboratory and research institute.
- Ion chromatography and titrator
- Pilot scale and industry production

Suitable Pump Head



Multichannel
Pump Head
(MC1-MC12)



Multichannel
Pump Head
(DG1-DG12)




Easy Load
Pump Head
(YZ1515x)



Easy Load
Pump Head
(YZ2515x)

LabN Series External Control Schematic Diagram

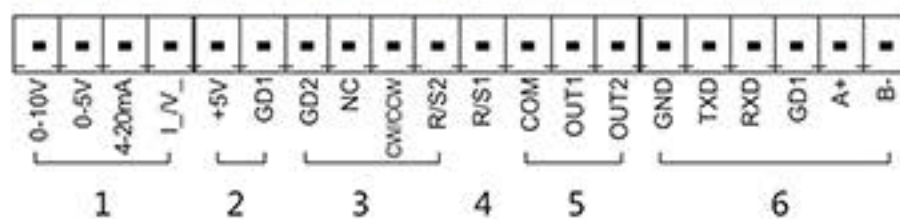

LabN Series Membrane Functions

**SHENCHEN****Standard Peristaltic Pump****Standard Peristaltic Pump**

BT100N,BT300N,BT600N

**Working
quiet****Features**

- Plastic coated metal housing, compact structure.
- Has timing function, time range 0.5s-999s, can be used for simple dispense.
- **RS232、RS485 Communication Interface.**

Support Shenchen communication protocol or standard Modbus communication protocol(RTU mode).



Model Number ■ BT100N ■ BT300N ■ BT600N

LabN Series/BT-N Series Technical Specifications

Flow rate range	LabN1: 0.000166~570 mL/min	Control method	Membrane keypad
	LabN3: 0.000166~1330 mL/min	External speed control signal	0~5V,0~10V,4~20mA for option
	LabN6: 0.000166~2280 mL/min	Start/stop,	Passive switch signal, such as foot pedal
	BT100N: 0.000829~570 mL/min	reversing signal	Active switch signal: 5V,12V,24V for option
	BT300N: 0.000829~1330 mL/min	Communication interface	Rs232, RS485 support Modbus protocol(RTU mode)
	BT600N: 0.000829~2280 mL/min	Output interface	LabN Series: Output motor working status (Open-Collector output)
Speed resolution	LabN Series :0.1rpm;	Power supply	AC 220V ± 10% 50Hz/60Hz (standard)
	BT-N Series: 0~100rpm,0.1rpm; 100~600rpm, 1rpm		AC 110V ± 10% 50Hz/60Hz (optional)
Back suction angle	LabN Series: 0~360°	Drive dimension	LabN Series: 261.4 × 157.3 × 236.9mm
Testing time range	LabN Series: 0.1 s~9999 h		BT-N Series: 183 × 131 × 194mm
	BT-N Series: 0.5 s~999 s	Drive weight	LabN Series: 4.40 kg BT-N Series:4.20 kg
Outlet pressure	0.86~1.0 mm wall thickness tubing: 0.1Mpa;	Power consumption	<50W
	1.6~2.4mm wall thickness tubing:0.1~0.27Mpa	Condition temperature	0~40℃
Display	LabN Series: 3.2" high definition LCD screen	Relative humidity	< 80%
	BT-N Series: LED Display	IP rate	IP31

LabN Series/BT-N Series Product Composition and Flow Rate Range

Standard peristaltic pump		Pump head & Flow rates (mL/min)			
		YZ1515x	YZ2515x	MC1~MC12(10)	MC1~MC12(6)
Drive & speed	Tubing	13', 14', 19', 16' 25', 17', 18'	15', 24'	Inner diameter: 0.13-3.17mm Wall thickness: 0.86-1mm	
LabN1	0.1-150rpm	0.007~570	0.17~435	0.000166-49(suggest working rotating speed ≤ 150rpm)	0.000185-65(suggest working rotating speed ≤ 150rpm)
LabN3	0.1-350rpm	0.007~1330	0.17~1015		
LabN6	0.1-600rpm	0.007~2280	0.17~1740		
BT100N	0.5-150rpm	0.035~570	0.85~435	0.000829-49(suggest working rotating speed ≤ 150rpm)	0.000925-65(suggest working rotating speed ≤ 150rpm)
BT300N	0.5-350rpm	0.035~1330	0.85~1015		
BT600N	0.5-600rpm	0.035~2280	0.85~1740		



Standard Peristaltic Pump



Features

- Servo motor drive, powerful and maintenance-free.
- 3 digital LED display motor speed, rotary encoded switch control.
- Suitable for industrial sites, transfer liquid with large flow and high precision.

Model Number

- N6-3L/DZ25-3L
- N6-6L/DZ25-6L
- N6-12L/YZ35

Technical Specifications

Flow rate range	N6-3L: 0.211~3600 mL/min N6-6L: 0.3~6000 mL/min N6-12L: 0.00069~12 L/min	Communication interface	RS232(standard),Shenchen-N232 Communication protocol; RS485(optional), Modbus protocol(RTU mode)
Speed range	0.1~600 rpm	Power supply	AC 220V $\pm$ 10% 50Hz/60Hz (standard) AC 110V $\pm$ 10% 50Hz/60Hz (optional)
Speed resolution	0~100rpm, 0.1rpm; 100~600rpm, 1rpm.	Drive dimension (L x W x H)	N6-3L: 223 x 152 x 168mm N6-6L: 283 x 192 x 218mm N6-12L: 302 x 222 x 275mm
Flow rate accuracy	$\pm$ 0.5%	Drive weight	N6-3L: 5.02kg; N6-6L: 7.85kg N6-12L: 13.14kg
Outlet pressure	0.3Mpa	Power consumption	N6-3L: <80W ; N6-6L: <180W N6-12L: <300W
Motor type	Servo motor	Condition temperature	0~40°C
Display	3 digital LED	Relative humidity	< 80%
Control method	Light touch membrane keypad	IP rate	IP31
Speed control	Rotary encoded switch		
Start/stop, direction signal	Passive switch signal, such as foot pedal Active switch signal: 24V default		
External speed control signal	0~5V(standard), 0~10V, 4~20mA(optional)		

Product Composition and Flow Rate Range

Drive	Motor Type	Pump Head	Tubing Size	Speed Range	Flow Rate Range(mL/min)
N6-3L	Servo Motor	DZ25-3L	15', 24', 35', 36'	0.1~600 (rpm)	0.211~3600
N6-6L		DZ25-6L	15', 24', 35', 36'		0.3~6000
N6-12L		YZ35	26', 73', 82'		0.69~12000

BAODING PERISTALTIC PUMP ENGINEERING TECHNOLOGY RESEARCH CENTER

New Arrivals

LabM series

ABS engineering
plastic housing



RESEARCH CENTER

BAODING SHENCHEN PRECISION PUMP CO.,LTD

Baoding Shencheng Precision Pump Co.,Ltd
Baoding Peristaltic Pump Engineering Technology Research Center



Add: No.103,Building 2, ZhiDian Industrial Park,FuXing East Road 999, Baoding,China.
Website: www.good-pump.com Tel : 0086-312-6780680 Fax : 0086-312-6780636



Basic Peristaltic Pump

LabM1, LabM3, LabM6



Features

- 3 digital LED display motor speed, membrane keypad control.
- Has timing function, can be used for simple dispense
- Has both RS232 and RS485 communication mode.

Model Number

- LabM1
- LabM3
- LabM6

Suitable Pump Head



Multichannel
Pump Head
(MC1-MC12)



Multichannel
Pump Head
(DG1-DG12)



Easy Load
Pump Head
(YZ1515x)



Easy Load
Pump Head
(YZ2515x)

Typical Application

- University laboratory.
- Supporting ion chromatography and titrator.



LabM Series Membrane keypad Schematic Diagram



SHENCHEN Basic Peristaltic Pump

Basic Peristaltic Pump

BT100M, BT300M, BT600M



Easy
operation

Features

- 3 digital LED display motor speed.
- Toggle switch control start/stop and direction.
- Rotary encoded switch control speed, easy operation.
- Economical and durable, suitable for different industrial application.

Model Number

- BT100M
- BT300M
- BT600M

Typical Application

- Add cement grinding aids
- Support industrial production line

LabM Series/BT-M Series Technical Specifications

Flow rate range	LabM1: 0.000829~570 mL/min LabM3: 0.000829~1330 mL/min LabM6: 0.000829~2280 mL/min BT100M: 0.00166~570 mL/min BT300M: 0.00166~1330 mL/min BT600M: 0.00166~2280 mL/min	Start/stop, direction signal	LabM Series: Passive switch signal, such as foot pedal switch; Active switch signal: 24V default BT-M Series: TTL level signal
Speed resolution	0.1rpm when the speed is 0~100rpm, 1rpm when the speed is 100~600rpm.	Communication interface	LabM Series: RS232(standard), Shenchen communication protocol; RS485(optional), Modbus protocol(RTU mode)
Testing time range	LabM Series: 0.5s~999s	Power supply	AC 220V $\pm$ 10% 50Hz/60Hz (standard) AC 110V $\pm$ 10% 50Hz/60Hz (optional)
Outlet pressure	0.1Mpa (0.86~1.0mm wall thickness tubing) 0.1~0.27Mpa (1.6~2.4mm wall thickness tubing)	Drive dimension	LabM Series: 261.4 $\times$ 157.3 $\times$ 236.9mm BT100M: 234 $\times$ 144 $\times$ 140mm BT300M, BT600M: 254 $\times$ 185 $\times$ 149mm
Back suction angle	BT-M Series: 0~360° (optional)	Drive weight	LabM Series: 4.40 kg BT100M: 3.6kg BT300M, BT600M: 4.64kg
Display	3 digital LED	Power consumption	<50W
Control method	LabM Series: Membrane keypad BT-M Series: Toggle switch+Rotary encoded switch	Condition temperature	0~40°C
External speed control signal	LabM Series: 0~5V, 0~10V, 4~20mA for option BT-M Series: 0~5V(standard), 0~10V, 4~20mA(option)	Relative humidity	<80%
		IP rate	IP31

LabM Series/BT-M Series Product Composition and Flow Rate Range

Basic peristaltic pump		Pump head & Flow rates (mL/min)			
		YZ1515x	YZ2515x	MC1-MC12(10)	MC1-MC12(6)
Drive & speed	Tubing	13', 14', 19', 16' 25', 17', 18'	15', 24'	Inner diameter: 0.13~3.17mm Wall thickness: 0.86~1mm	
LabM1	0.5~150rpm	0.035~570	0.85~435	0.000829~49(suggest working rotating speed $\leq$ 150rpm)	0.000925~65(suggest working rotating speed $\leq$ 150rpm)
LabM3	0.5~350rpm	0.035~1330	0.85~1015		
LabM6	0.5~600rpm	0.035~2280	0.85~1740		
BT100M	1~150rpm	0.07~570	1.7~435	0.00166~49(suggest working rotating speed $\leq$ 150rpm)	0.00185~65(suggest working rotating speed $\leq$ 150rpm)
BT300M	1~350rpm	0.07~1330	1.7~1015		
BT600M	1~600rpm	0.07~2280	1.7~1740		




Explosion Proof Peristaltic Pump EXP600/YZ35



Features

- Driven by explosion-proof motor, the power is strong.
- Stepless speed regulation, reversible.
- IP rate: IP55
- Suitable for inflammable and explosive industrial sites.

Installation Sizes Drawing (Unit: mm)

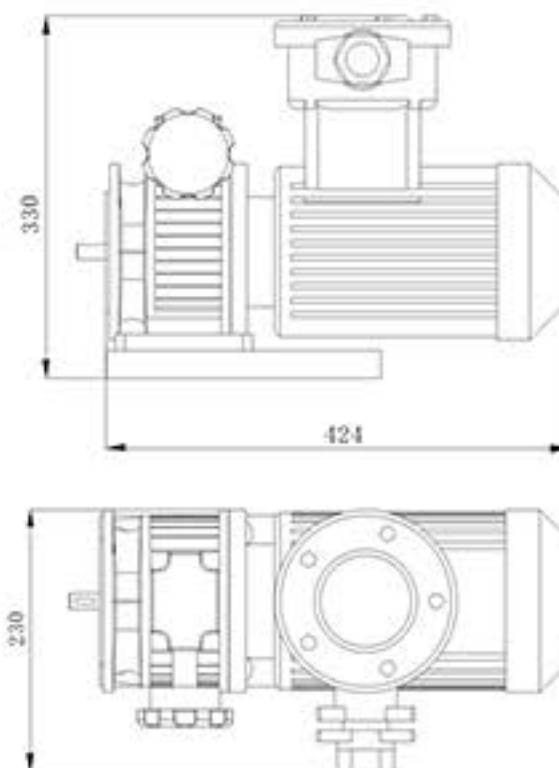

EXP600/YZ1515x



Technical Specifications

Speed range	180-600 rpm
Speed control	Manual stepless speed regulation
Voltage	AC 3 phase 380V
Condition temperature	0-40°C
Drive dimension	424 x 230 x 330 (mm) (L x W x H)
Drive weight	30 kg
Relative humidity	<80%
IP rate	IP55

Product Composition and Flow Rate Range

Drive	Speed	Pump Head	Tubing Size	Flow rates (mL/min)
EXP600	180-600rpm	YZ1515x	13', 14', 19', 16', 25', 17', 18'	12.6-2280
		YZ2515x	15', 24'	306-1740
		YZ35	26', 73', 82'	1242-12000
		DZ25-3L	15', 24', 35', 36'	379.8-3600
		DZ25-6L	15', 24', 35', 36'	540-6000
		SN15	14', 16'	43.2-528
		SN25	24'	450-1500



SHENCHEN Explosion Proof Peristaltic Pump

Explosion Proof Peristaltic Pump

QD600

NEW



Features

- Pneumatic motor driver, gas driving, explosion-proof and safe.
- Stepless speed regulation, manual control flow valve.
- Start/stop, direction, overload protection functions.
- Barometer, convenient for users to check the pressure value.
- Suitable for inflammable, explosive industrial sites.

Suitable Pump Head



Multichannel
Pump Head
(MC1-MC12)



Multichannel
Pump Head
(DG1-DG12)




Easy Load
Pump Head
(YZ1515x)




Easy Load
Pump Head
(YZ2515x)

Installation Sizes Drawing (Unit: mm)

Working status dimension drawing

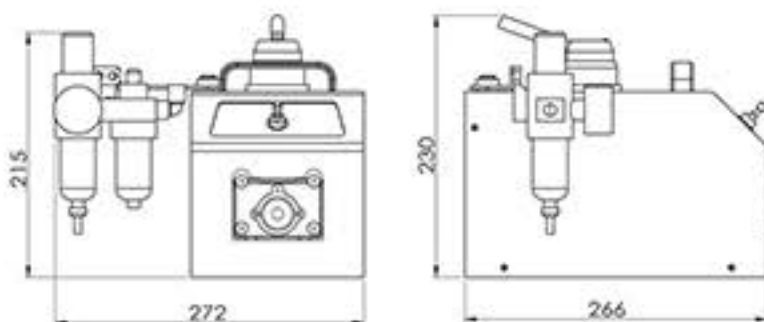

Non-working status dimension drawing

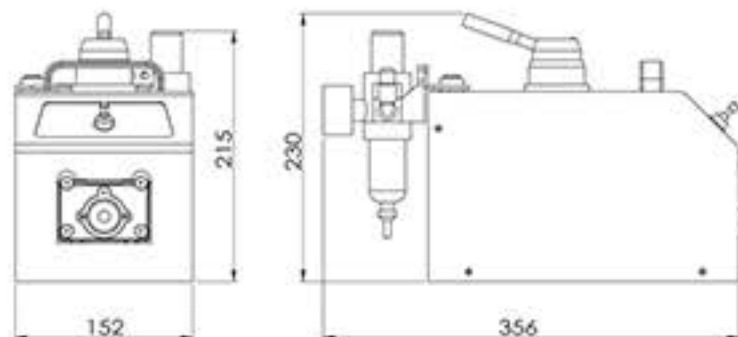

Technical Specification

Speed range	60~600 rpm
Display	Display air pressure
Speed control	stepless speed regulation (By adjusting the flow valve)
Control function	Start/stop, reversing (gas reversing valve control)
Working air pressure	0.1~0.4Mpa
Gas consumption	4L/Sec
Drive dimension (L x W x H)	Working: 266x 272 x 230mm Non-working: 356x 152 x 230mm
Power consumption	< 80W
Condition temperature	0~40°C
Relative humidity	< 80%
IP rate	IP31
Drive weight	4.10 kg

Product Composition and Flow Rate Range

Peristaltic Pump		Pump head & Flow rates (mL/min)			
		YZ1515x	YZ2515x	MC1-MC12(10)	MC1-MC12(6)
Drive & speed	Tubing	13', 14', 19', 16' 25', 17', 18'	15', 24'	Inner diameter: 0.13~3.17mm Wall thickness: 0.86~1mm	
QD600	60-600rpm	5~2280	100~1740	2.8~49(suggest working rotating speed ≤ 150rpm)	3.2~65(suggest working rotating speed ≤ 150rpm)



Medical Dialysis Pump UC15




Medical
Dialysis Pump
NEW

Product Introduction


UC15 spring pump head is widely used for medical dialysis equipment. Transparent protective cover, can observe the working situation anytime; uncover signal output, stop working when uncover to ensure safety; Elastic roller assembly, running smoothly; manual lever help install the tube, and can finish liquid circulation when power down unexpectedly; support stepper motor or brushless DC motor.

Installation Sizes Drawing (Unit: mm)

Support brushless DC motor

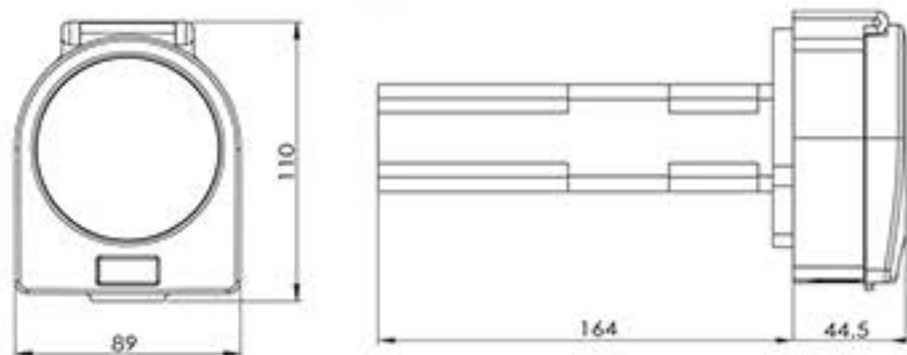

Uncover output signal description

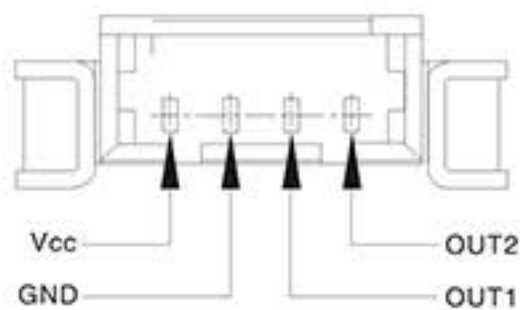

Motor data

Speed range	0.1-200 rpm
Speed resolution	0.1 rpm
Motor type	Stepper motor
Output signal	Uncover signal (Active signal DC12V)



1、Press the pipe clamp, put the tube into the U shaped groove.



2、Put the tube between the two guide posts. Hold the spanner, rotate the wheel body.



3、Press the pipe clamp, put the tube into the U shaped groove.



4、Installation is finished.


Technical Specifications

Model	Housing material	Tubing		Speed range(rpm)	Flow rates (mL/min)	Weight(kg)
		Size	ID x Wall(mm)			
UC15	Aluminum alloy /PC	19#	2.4 x 1.6	0.1-200	0.08-183	0.5
		16#	3.1 x 1.6		0.15-320	
		25#	4.8 x 1.6		1.1-699	
		17#	6.4 x 1.6		2-1270	
		18#	7.9 x 1.6		2.8-1719	
		Medical dialysis tube		0.1-60	0.86-500	



SHENCHEN Compact Pump Head

Compact Pump Head MiniPump Series

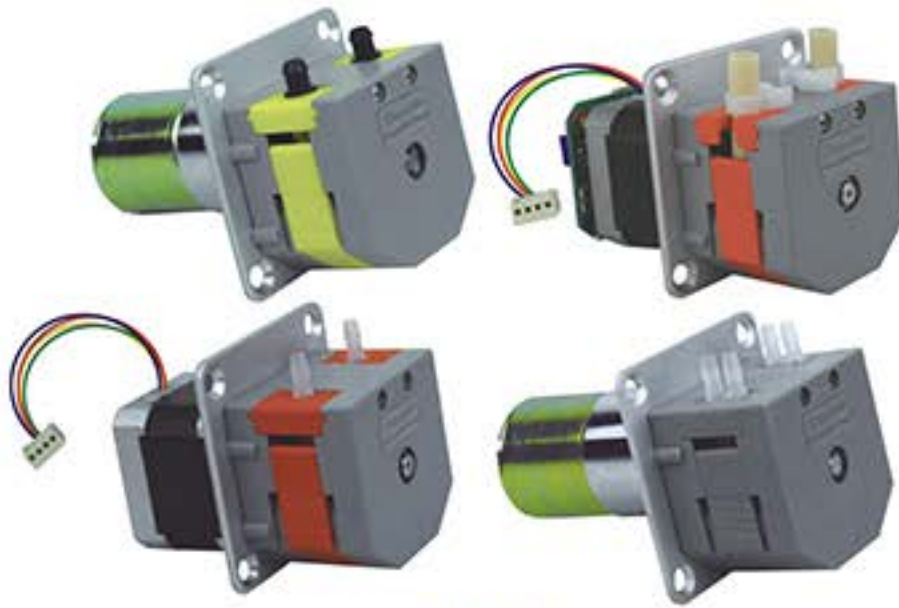

Patent No.: ZL2011 3 0002460.5
Patent No.: ZL2011 2 0015673.6



Tubing Installation Procedure

A. Make the two tubing clamps align the slide way of pump head two sides.



B. Ensure that the tubing is in the middle of rollers, and push the clamps into the pump head.



C. Make the compression block align the upper slide way of pump head.



D. Press the compression block into the pump head, and ensure it fixed.



Product Introduction

Minipump series products which adopt compact and embedded design are more suitable for supporting a variety of instruments and equipments. Single and dual channel are optional. It accepts 8 kinds of tubing. **Flow rates are 0.0024-190 mL/min.** PharMed BPT long-life tubing from French Saint-Gobain is installed in internal part of pump head. Applicable motors for Minipump are stepper motor, DC motor, AC motor and other modes to drive with wide applications.

Typical Application


- Discharge condensate water in the flue gases survey meter.
- Support biochemistry analytical instruments.

Features


- Attractive appearance, compact structure.
- Easy to operate and economical cost.
- Accept many kinds of motors to drive.
- Supply several colors of appearance, ideal for supporting analytical instruments.



Installation Sizes Drawing (Unit: mm)

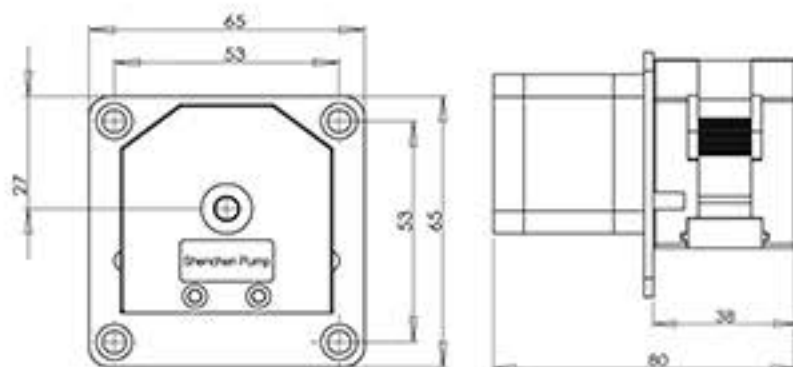

MiniPump-Stepper motor drive

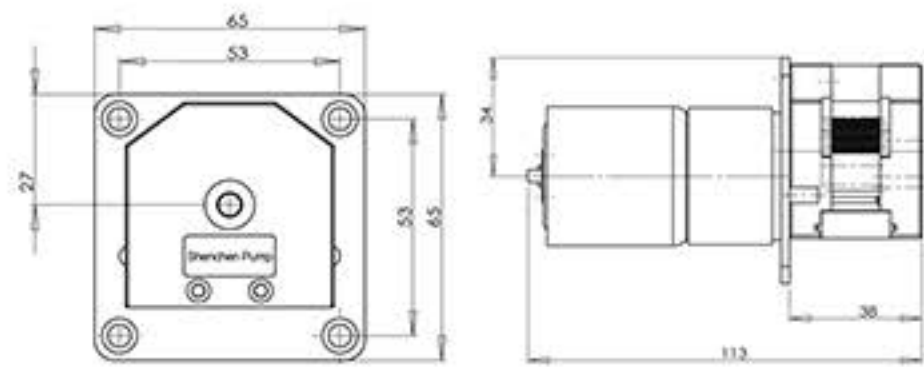



MiniPump-DC motor drive



Compact Pump Head MiniPump-N Series



Installation Sizes Drawing (Unit: mm)

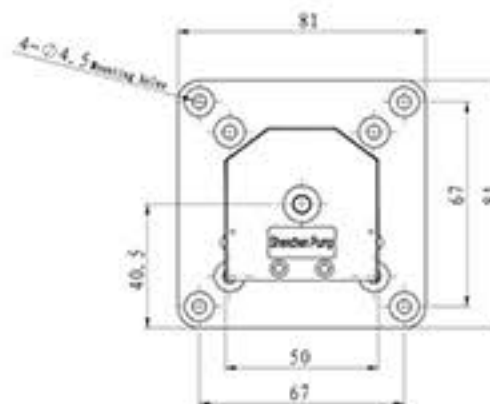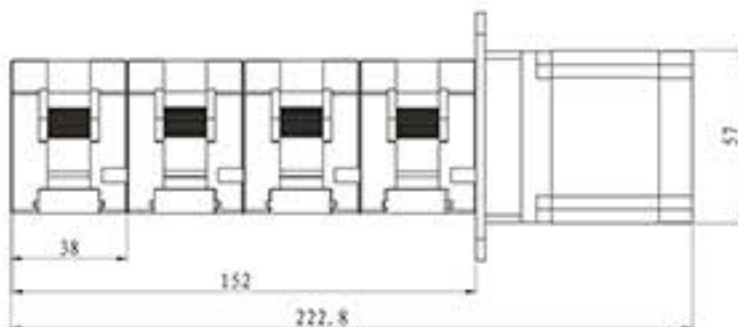

More details:

MiniPump-N series products can adopt single and dual channel pump head, suitable for 8 different sizes tube, use stepper motor drive, can combine a variety of colors.

MiniPump Series and MiniPump-N Series Technical Specifications

Rotation speed (rpm)	Dual channel Flow rates (mL/min) and tubing (1 mm wall thickness)			
	ID 1 mm	ID 2 mm	ID 2.5 mm	ID 3 mm
5	0.25	0.90	1.28	1.78
10	0.49	1.84	2.59	3.52
15	0.74	2.76	3.85	5.30
25	1.25	4.57	6.43	8.80
40	2.00	7.31	10.25	14.02
50	2.52	9.14	12.78	17.47
75	3.76	13.64	19.04	26.12
100	5.03	18.17	25.42	35.34
200	10.12	36.45	51.27	71.44
250	12.63	45.75	64.34	88.35
300	15.01	54.63	76.84	108.39

Rotation speed (rpm)	Single channel Flow rates (mL/min) and tubing (1.6 mm wall thickness)			
	13# Tubing ID 0.8 mm	14# Tubing ID 1.6 mm	16# Tubing ID 3.1 mm	25# Tubing ID 4.8 mm
5	0.12	0.56	1.97	3.26
10	0.26	1.15	3.90	6.41
15	0.39	1.74	5.80	9.52
25	0.66	2.88	9.59	15.74
40	1.08	4.59	15.29	25.15
50	1.34	5.69	18.96	31.14
75	2.00	8.54	28.36	46.55
100	2.73	11.37	37.66	61.98
200	5.46	22.76	76.02	124.14
250	6.81	28.40	95.08	153.63
300	8.28	33.88	114.31	190.00



SHENCHEN Multichannel Pump Head

Multichannel Pump Head DG Series




Product Introduction

The first generation multichannel pump head is used for multichannel transfer micro liquid. Occlusion can be adjusted. Widely used in supporting analytical instruments.

Model Number

- DG-1 ■ DG-2 ■ DG-4
- DG-6 ■ DG-8 ■ DG-12

Installation Sizes Drawing (Unit: mm)

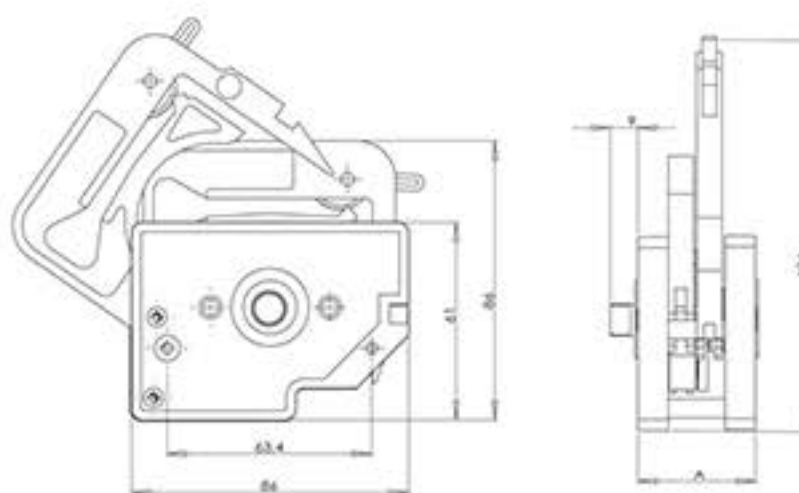

Typical Application

- Widely used in supporting analytical instruments and laboratories.



DG-4, DG-6, DG-8, DG-12 use unique spanner to open the cartridge

Technical Specifications

Pump head		DG-1	DG-2	DG-4	DG-6	DG-8	DG-12
Channel number		1	2	4	6	8	12
Tubing		Inside Diameter: 0.13–2.5mm, wall thickness: 0.5–1mm					
Flow rates (≤100rpm)	6 rollers	0.00008–27 mL/min(per channel)					
	10 rollers	0.000067–21.6 mL/min (per channel)					
Cartridge material		POM					
Rollers material		304 stainless steel					
Weight (kg)	6 rollers	0.21	0.27	0.40	0.54	0.67	0.95
	10 rollers	0.20	0.26	0.39	0.51	0.63	0.88



Multichannel Pump Head-MC Series POM Material & PVDF Material

Appearance design patent No.: ZL2011 3 0101376.9
The utility model patent No.: ZL2011 2 0189392.2




To open the cartridge from upside will be more convenient.

Product Introduction

MC series new generation multichannel pump head with unique cam structure, make it much easier to operate. It is designed for small flow rates and multichannel transfer. Comparing with traditional DG pump head, it has particular cam design for micro adjustment occlusion, and also for opening the cartridge from the top. 304 stainless steel rollers assembly include two structures of 6 rollers and 10 rollers. Install silicone tubing and TYGON series imported tubing. It is widely used in supporting analytical instruments.

Installation Sizes Drawing (Unit: mm)

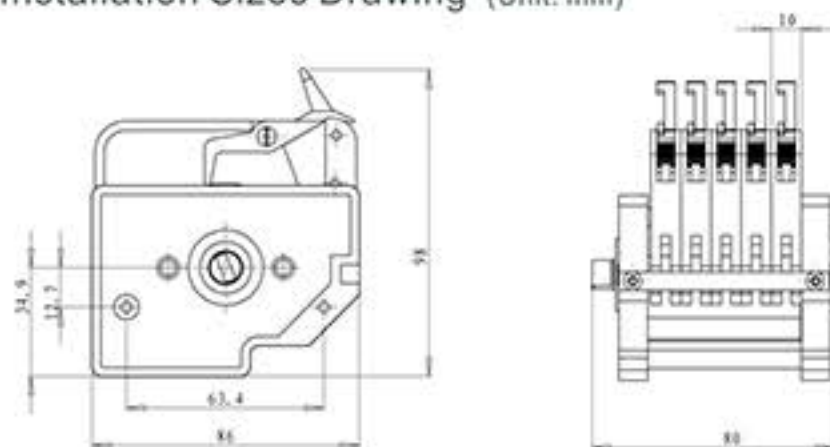

Typical Application


- Support blood sugar analytical device.
- Support ammonian analytical instrument.

Note: Increasing or reducing one channel, the pump head size need to increase or reduce the thickness of one cartridge.



Closed design

Installation Sizes Drawing (Unit: mm)

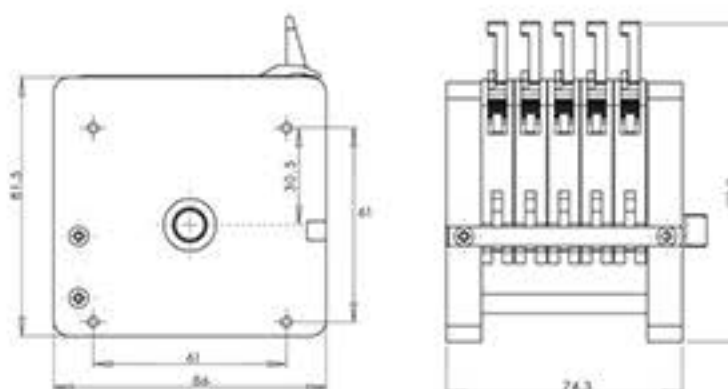

Note: Increasing or reducing one channel, the pump head size need to increase or reduce the thickness of one cartridge.

MC Pump VS DG Pump



1. The traditional DG series pump head need to open the cartridge from the bottom. It is very inconvenient to use when the space of the equipment is compact. The MC series pump head just need to open the cartridge from the top, easy to operate.


2. The gear on the MC pump head is smaller than DG, this improved the adjusting resolution.



DG Series



SHENCHEN Multichannel Pump Head

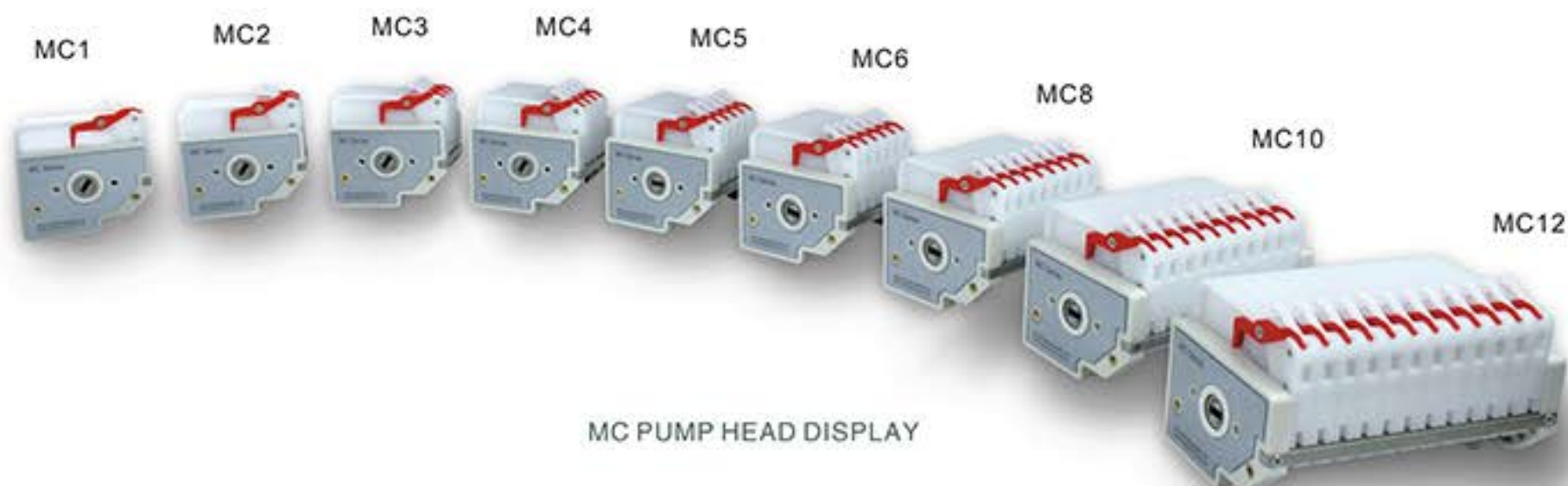

Tubing Installation Procedure

A. Turn the spanner clockwise to take it off from the lever.



B. After that, the cartridge can be taken out lightly.



C. Install the tubing and confirm the working position between two adjacent clamps.



D. Embed the cartridge and press it lightly to reach the working position.



E. Turn the spanner counter-clockwise to make it reach the lever.



F. Multichannel pump heads and tubing drawing




MC Technical Specifications

Channel number	1	2	3	4	5	6	8	10	12
Cartridge material	POM/PVDF								
Rollers material	304/316 stainless steel material								
Rollers number	6 rollers, 10 rollers,								

Tubing	Speed	Flow rates of pump head with 10 rollers (mL/min)	Flow rates of pump head with 6 rollers (mL/min)	Tubing maximum pressure (Mpa)	
1×1	0.1~150rpm	0.0046~7.03	0.0053~8.13	Continuous	0.1
2×1		0.016~23.07	0.02~33.08	Intermittent	0.1
2.4×0.8		0.021~31.26	0.027~44.62		
3×1		0.032~49.07	0.043~65.00		



Multichannel Pump Head MD3

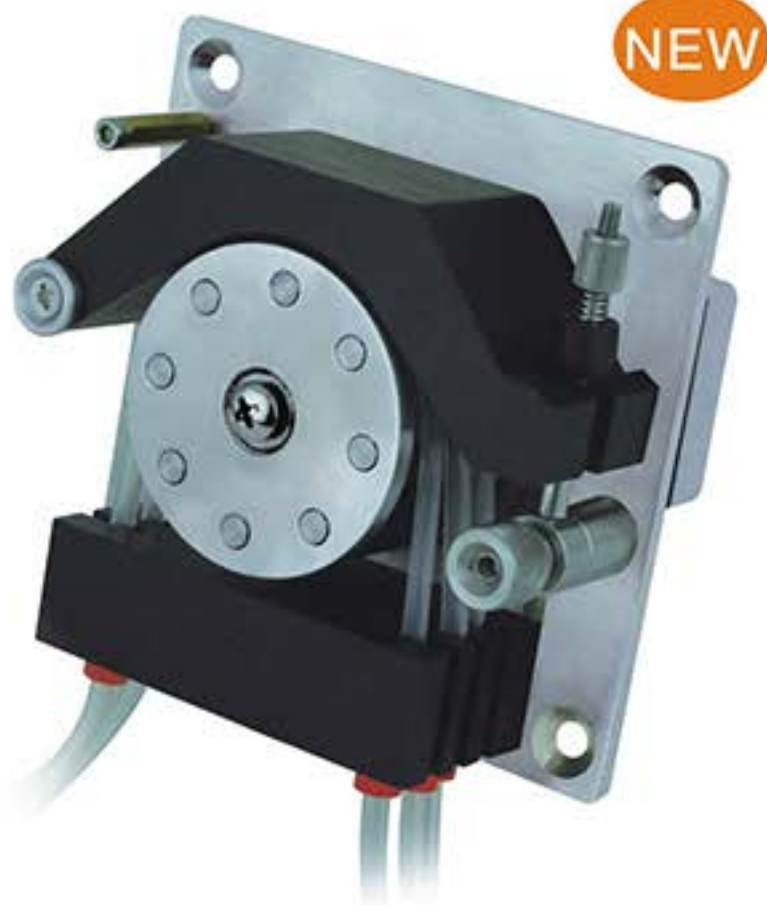

Model Number ■ MD3


Typical Application

- Support analytical instruments

Product Introduction

Suitable for liquid transferring of small flow rate and multiple channels. Complete compression block design, can improve the accuracy among channels effectively; with the function of micro adjustment occlusion, to micro adjust the flow rate and outlet pressure. Eight rollers design, reduce the liquid pulsation; embedded design, very suitable for equipment supporting.

Installation Sizes Drawing (Unit: mm)

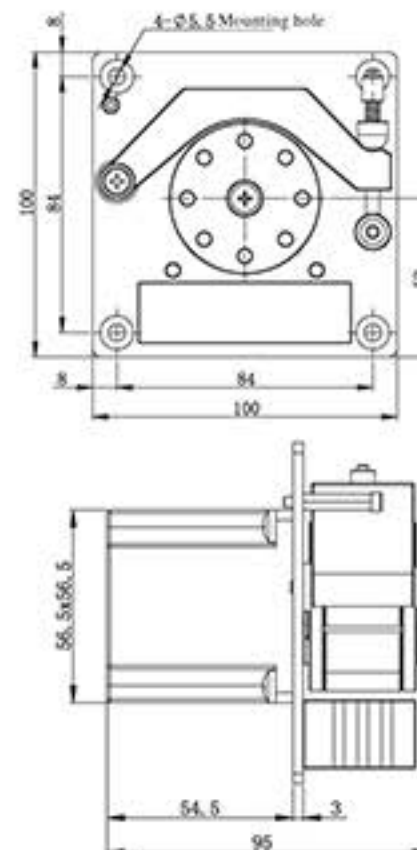

Tubing Installation Procedure



1. Turn the spanner clockwise.



2. Install tubing one by one, from inside to outside.



3. Tubing installation finished.



4. Turn spanner anticlockwise, locking the upper block, micro adjust the occlusion.

Technical Specifications

Pump head	Channel number	Speed	Tubing	Flow rates (0.1-300rpm)	Rollers number	Rollers material	Cartridge material	Tubing maximum pressure (Mpa)	
MD3	3	0.1-300 (rpm)	1 × 1	0.007-20.95 mL/min (each channel)	8 rollers	304 stainless steel	POM	Intermittent 0.1	Continuous 0.1
			2 × 1	0.0244-74mL/min (each channel)					
			2.79 × 0.86	0.0356-107mL/min (each channel)					




SHENCHEN Multichannel Pump Head

Multichannel Pump Head

MT12



Installation Sizes Drawing (Unit: mm)

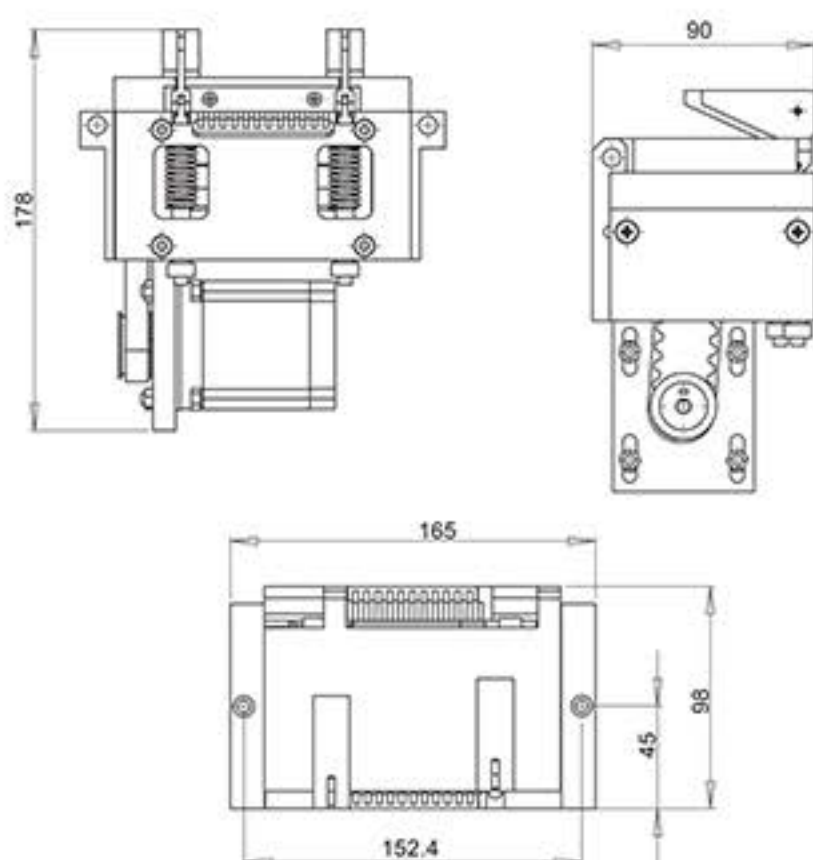

Product Introduction

- The unique complete compression solve the big error problem between channels from traditional multichannel pump.
- Accept many kinds of motors drive, such as stepper motor, DC motor, AC motor and servo motor.

Model Number

- MT12

Typical Application

- Support the chemistry analytical instruments

Technical Specifications

Model No.	Housing material	Rollers material	Channel number	Rollers number	Speed (rpm)	Tubing and flow rates per channel (mL/min)				Tubing maximum pressure (Mpa)	
						0.51 × 0.9	1.02 × 0.86	2.06 × 0.86	2.54 × 0.86	Intermittent	Continuous
MT12	Aluminum alloy	304/316 stainless steel	12 channels	10 rollers	0.1-50	0.0016-0.84	0.005-2.5	0.015-7.5	0.020-10.5	0.2	0.15



Easy Load Pump Head

YZ1515x, YZ2515x

PSF




PPS



Product Introduction

Change the tubing easily and rapidly; automatic tubing retention; accept 9 sizes tubing to meet different flow rates requests. Two kinds of material housing for customers' selection. PSF (polysulphone material): white, exquisite appearance, good rigidity, high precision, usually used in research laboratory and supporting analytical device. PPS (polyphenylene sulfide material): black, resist organic solvent corrosion, high precision and rigidity, suitable for industry or other bad working conditions. Rollers assembly with 304 stainless steel material resists corrosion with long lifetime.

Installation Sizes Drawing (Unit: mm)

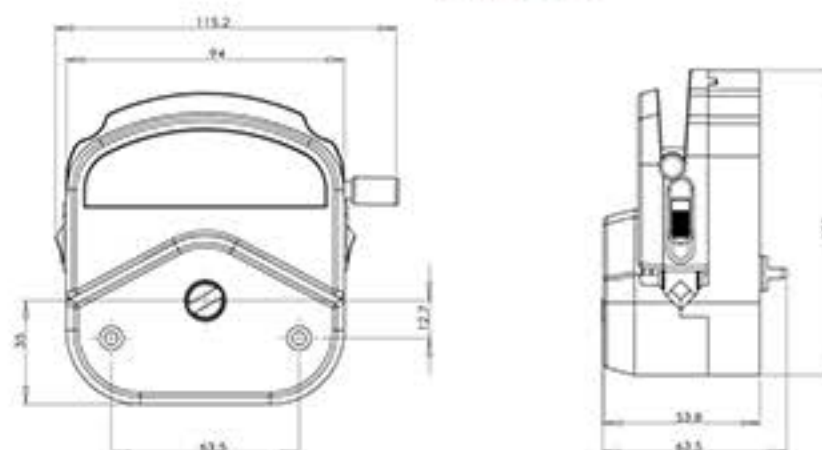

Mounting several pump heads can increase the flow rates.



Model Number ■ YZ1515x ■ YZ2515x

Typical Application

- Support COD on-line monitor.
- Support anorectal therapeutic instrument.
- Support biofermentation cylinder.



3 rollers

6 rollers

YZ Series pump head has wide application, and is usually used together with many kinds of analytical instruments. According to the different technical requests, it can adopt different motors to drive, such as stepper motor, DC motor, synchronous motor, AC gear motor, etc. Stack up several pump heads can increase the flow rates.

Pump head	Tubing	ID×Wall thickness(mm)	mL / r	Flow rates(mL/min) (1-600rpm)	Tubing Max.pressure(Mpa)		Material/Weight(kg)	
					Intermittent	Continuous	PSF	PPS
YZ1515x	13"	0.8×1.6	0.07	0.007~42	0.27	0.17	0.40 (3 rollers)	0.46 (3 rollers)
	14"	1.6×1.6	0.27	0.027~162				
	19"	2.4×1.6	0.55	0.055~330				
	16"	3.1×1.6	0.82	0.082~492				
	25"	4.8×1.6	1.7	0.17~1020	0.24	0.14	0.44 (6 rollers)	0.50 (6 rollers)
	17"	6.4×1.6	2.9	0.29~1740	0.14	0.10		
	18"	7.9×1.6	3.8	0.38~2280	0.10	0.07		
YZ2515x	15"	4.8×2.4	1.7	0.17~1020	0.27	0.17		
	24"	6.4×2.4	2.9	0.29~1740				



SHENCHEN Easy Load Pump Head

Easy Load Pump Head

YZ35-PPS, YZ35(Aluminum Alloy)



PPS



Aluminum Alloy



3 rollers

Product Introduction


YZ35 pump head with aluminum alloy or PPS housing material. 304 stainless steel rollers assembly achieve high precision transferring liquid. Adaptive tubing cartridge structure makes it more easy to load the tubing.

Model Number


- YZ35-PPS
- YZ35(Aluminum Alloy)

Installation Sizes Drawing (Unit: mm)

YZ35-PPS

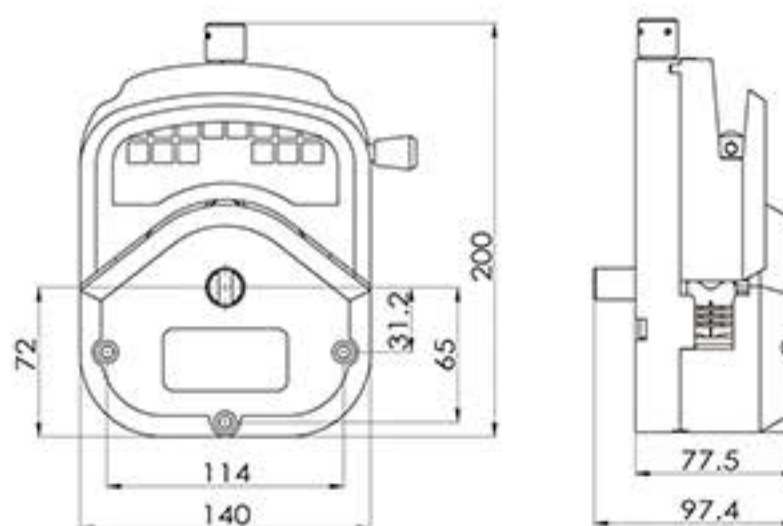

YZ35-Aluminum Alloy

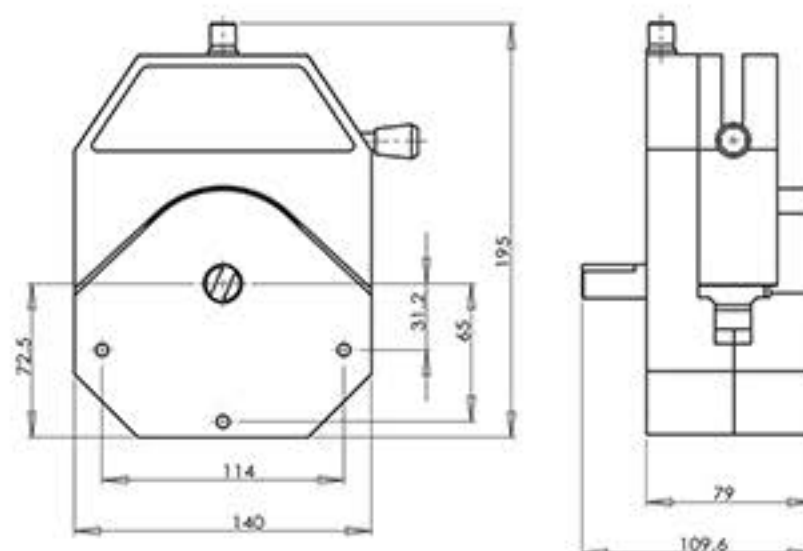

Technical Specifications

Pump head	Tubing	ID × Wall thickness(mm)	mL / r	Flow rates(mL/min) (1-600rpm)	Tubing Max.pressure(Mpa)		Material/Weight(kg)	
					Intermittent	Continuous	Aluminum Alloy	PPS
YZ35	26"	6.4 × 3.3	6.9	0.69~4200	0.27	0.20	4.36	1.50
	73"	9.6 × 3.3	12.3	1.23~7400				
	82"	12.7 × 3.3	20	2~12000	0.14	0.10		



Fast Load Pump Head

DZ-3L, DZ25-6L (Aluminum Alloy)




Product Introduction

Easy operation, change tube rapidly. The house material is aluminum alloy or PPS (polyphenylene sulfite). 304 stainless steel rollers. DZ25 pump head use 2.4mm wall thickness tube, the flow rate range is 0.2-6000mL/min. This pump head is suitable for high viscous liquid which include granule and floc. It also suitable for high pressure request. New cartridge structure design, can use the special connector, connect low cost tubing outside the pump head to save cost; also can use flexible tube card, meet hygiene requirement.

Model Number

■ DZ25-3L ■ DZ25-6L (Aluminum Alloy)

Installation Sizes Drawing (Unit: mm)

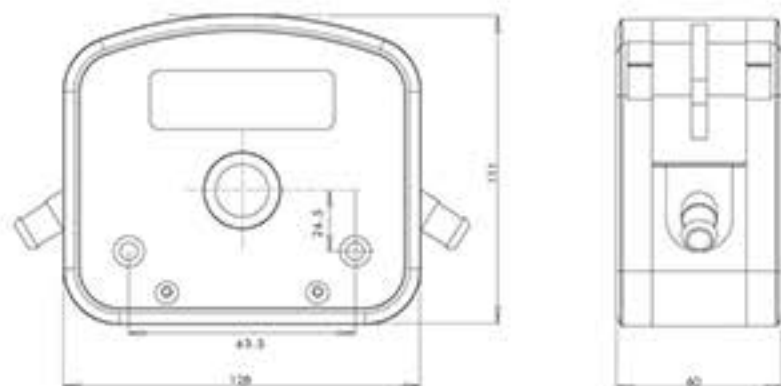

DZ25-6L Dimension



DZ25-3L Dimension

Tubing Installation Procedure


1. Lift both side levers, take off the upper block.



2. Put the tubing with cartridge or connector into pump housing.



3. Install the upper block, put down the levers to lock the block.

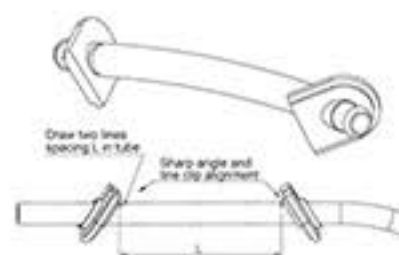

① Connector

The tubing length is 128mm between both connectors.

② Cartridge

DZ25-6L: the tubing length is 125mm between both cartridges.

DZ25-3L: the tubing length is 90mm between both cartridges.

Model No.	Housing Material	Cartridge/ Connector material	Tubing		Flow rates (mL/min) (0.1-600rpm)	Weight (kg)
			Tubing sizes	ID x Wall thickness		
DZ25-6L	Aluminum alloy	PP	15"	4.8 x 2.4 (mm)	0.3-1800	1.86
			24"	6.4 x 2.4 (mm)	0.55-3300	
			35"	7.9 x 2.4 (mm)	0.8-4800	
			36"	9.6 x 2.4 (mm)	1-6000	
DZ25-3L	PPS/PSF	PP	15"	4.8 x 2.4 (mm)	0.211-1264	1.16/0.5
			24"	6.4 x 2.4 (mm)	0.385-2310	
			35"	7.9 x 2.4 (mm)	0.508-3050	
			36"	9.6 x 2.4 (mm)	0.6-3600	



SHENCHEN Low Pulsation Pump Head

Low Pulsation Pump Head DY15,DY25

Aluminum Alloy



Product Introduction

The 'low pulsation' pump head is specialized for high precision filling. Through the phase difference between the two sets of rollers, make the pulse crest and trough complementary, reducing the fluid pulsation effectively, reach high precision liquid transferring.

Model Number

■ DY15 ■ DY25

Typical Application

- High precision dosing micro liquid.

Tubing Installation Procedure

1. Lift both side levers, take off the upper block.




2. Put the tubing with cartridge or connector into pump housing.



3. Install the upper block, put down the levers to lock the block.



Pulsation Reduction Schematic Diagram

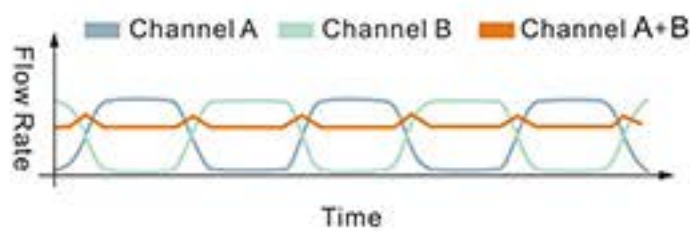


"Y" connector and tubing cartridge schematic diagram

The tubing length is 126mm between two cartridges



Installation Sizes Drawing (Unit: mm)

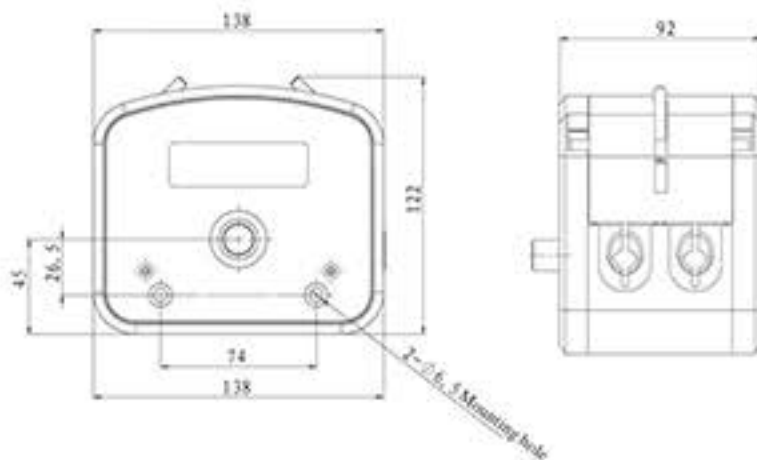

Technical Specifications

Model No.	Housing Material	Cartridge/Connector material	Tubing		Flow rates (mL/min) (0.1-350rpm)	Weight (kg)
			Tubing sizes	ID x Wall thickness		
DY15	Aluminum alloy	PP	13'	0.8 x 1.6(mm)	0.01~48	3.2
			14'	1.6 x 1.6(mm)	0.06~223	
			19'	2.4 x 1.6(mm)	0.13~448	
			16'	3.1 x 1.6(mm)	0.2~723	
			25'	4.8 x 1.6(mm)	0.47~1626	
			17'	6.4 x 1.6(mm)	0.64~2230	
			18'	7.9 x 1.6(mm)	0.95~3337	
DY25	Aluminum alloy	PP	15'	4.8 x 2.4(mm)	0.42~1480	3.2
			24'	6.4 x 2.4(mm)	0.76~2670	
			35'	7.9 x 2.4(mm)	0.1~360	
			36'	9.6 x 2.4(mm)	1.24~4340	



Standard Pump Head


SN25, SN15-14, SN15-16



Typical Application

- Use in the ice cream machine and coffee making machine
- Support COD on-line monitor

Installation Sizes Drawing (Unit: mm)

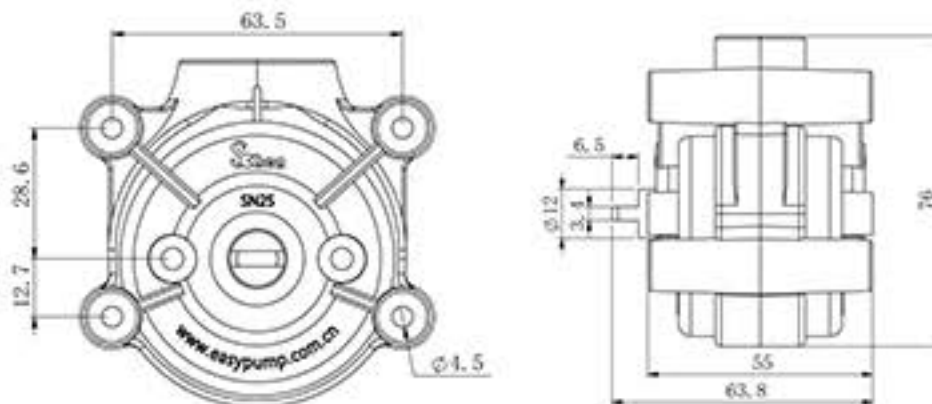

SN15-14



SN15-16

Product Introduction

Standard pump head with compact structure, exquisite appearance, fixed occlusion and stable running. Can stack several pump heads, and support various equipment.

Model Number

- SN25
- SN15-14
- SN15-16

Tubing Installation Procedure

A. Disassemble the pump head, and press the tubing on the rollers.



B. Make the tubing go around the rollers one cycle, and merge it on the exit.



C. Tighten the tubing aptly, and load the other half pump head body.



Technical Specifications

Model No.	Housing material	Rollers material	Tubing		Flow rates (mL/min) (0.1-600rpm)	Weight (kg)
			Tubing sizes	ID × Wall thickness		
SN15-14	PC	304 stainless steel	14"	1.6 × 1.6(mm)	0.024~144	0.20
SN15-16			16"	3.1 × 1.6(mm)	0.088~528	
SN25			24"	6.4 × 2.4(mm)	0.25~1500	



OEM Product
UD15



Tubing Installation Procedure



1. As below picture, turn the holding pin to unlocking position



2. Open the transparent cover



3. Install the tubing



4. Close the transparent cover, and turn the holding pin to unlocking position

Product Introduction

- PES material, transparent protecting cover, easily to observe the working situation of pump head and tubing at any time. Various of colors are available.
- There is water-proof groove in the pump head mounting side.
- Can use stepper motor, DC motor, AC motor and so on to drive.


Model Number

- UD15

Typical Application

- Support various of instrument

Installation Sizes Drawing (Unit: mm)

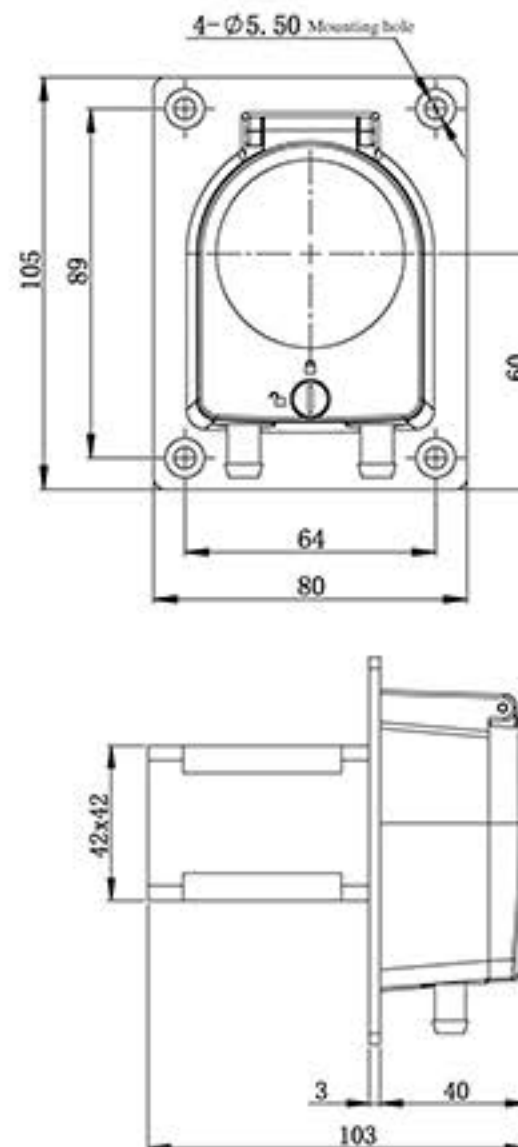

Technical Specifications

Model No.	Housing material		Tubing		Speed Range(rpm)	Flow rates (mL/min)	Weight (kg)
	Body	Protective cover	Tubing sizes	ID x Wall thickness			
UD15	PES	Transparent PC	16'	3.1 x 1.6(mm)	0.1~350	0.08~280	0.12
			25'	4.8 x 1.6(mm)		0.16~580	
			17'	6.4 x 1.6(mm)		0.30~1000	



OEM-B01




Outline dimensional drawing



Speed range	≤300 rpm
Max flow rate	190mL/min
Motor type	42 stepper motor
Motor phase number	2
Motor step angle	1.8
Motor phase voltage	2.2V
Motor phase current	1.0A

OEM-MZ Serie

Outline dimensional drawing

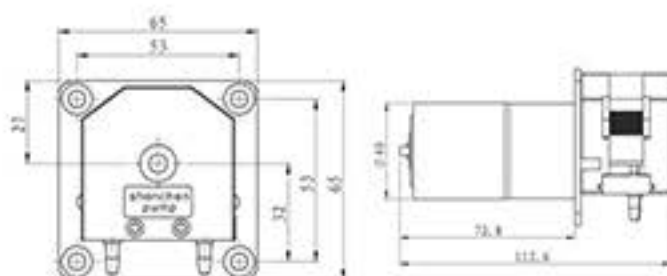


Model number	Rotation speed (rpm)	Power supply	Motor type
OEM-MZ01	170	DC 12V or DC 24V	Ø 37 DC Motor
OEM-MZ11	5		
OEM-MZ12	10		
OEM-MZ13	15		
OEM-MZ14	20		
OEM-MZ15	30		
OEM-MZ16	50		
OEM-MZ17	100		
OEM-MZ18	160		
OEM-MZ19	167		
OEM-MZ20	278		

OEM-MJ Serie

Outline dimensional drawing

NEW


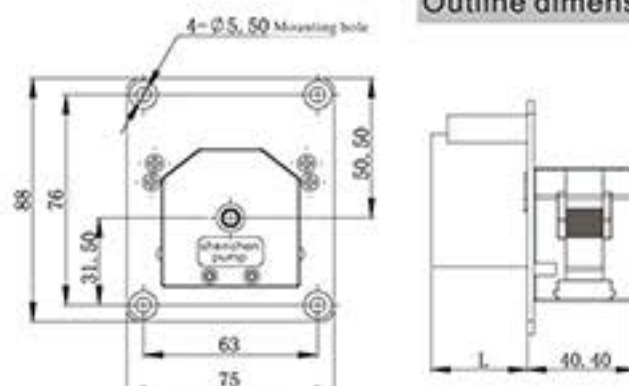


Model number	Rotation speed (rpm)	Direction of rotation	Power supply	Motor type
OEM-MJ101	2.5/3	Unidirectional	AC 110V or AC 220V	AC Motor
OEM-MJ102	4/5			
OEM-MJ103	5/6			
OEM-MJ201	4.2/5	Bidirectional		
OEM-MJ202	5/6			
OEM-MJ203	11.7/14			

Note: In the rotation of speed, "r" represents the rotation speed in 50/60HZ.

OEM-B126


Outline dimensional drawing

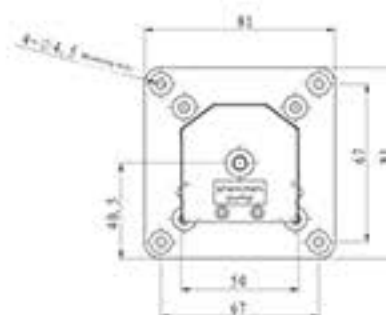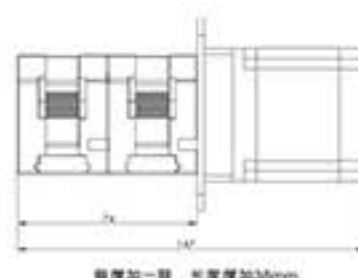

Speed range	≤300 rpm
Max flow rate	190mL/min
Pump head	MiniPump-N
Motor type	57 stepper motor
Motor phase number	2
Motor step angle	1.8
Motor phase voltage	3.75V
Motor phase current	3.0A

OEM-B013, OEM-B014

Outline dimensional drawing

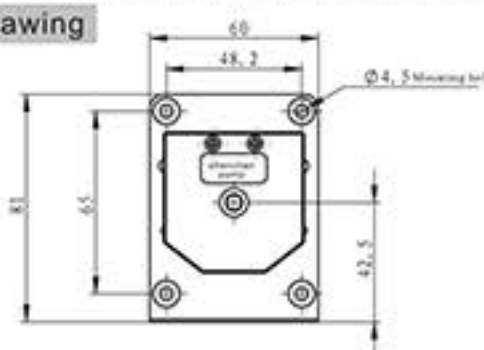

OEM-B013 Technical Specification

Speed range	0.1-300rpm
Flow rate range	0.0024-190ml/min
Communication interface	RS232, RS485 interface modbus protocol, control start/stop, direction and speed
Address setting	1-32
Motor type	42 stepper motor
Input voltage	DC 24V
Power	<30W

OEM-B014 Technical Specification


Speed range	0.1-300rpm
Flow rate range	0.12-190ml/min
Speed control	Different dial code to different speed 5-300rpm
External analog signal control	0-5V, 0-10V, 4-20mA 3 analog signal for option
External control	Start/stop, direction (external passive signal control)
Motor type	42 stepper motor
Input voltage	DC 24V
Power	<30W



OEM-B02、OEM-B03



Outline dimensional drawing

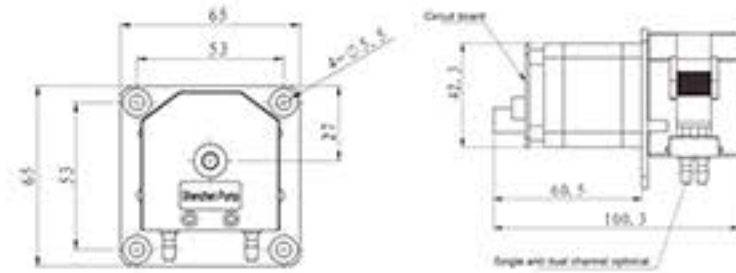

OEM-B02 Technical Specification

Speed range	0.1–300 rpm
Flow rate range	0.0024–190mL/min
Speed control	Preset speed or RS232 speed control
Subdivision settings	1/16
Motor type	42 stepper motor
Input voltage	DC 24V
Power	< 24W


OEM-B03 Technical Specification

Speed range	0.1–300 rpm
Flow rate range	0.0024–190mL/min
Speed control	External input pulse frequency speed control
Multiple segments	1/16(recommend), 1, 1/2, 1/4, 1/8, 1/32, 1/64, 1/128
Motor type	42 stepper motor
Input voltage	DC 24V
Power	< 24W

OEM-B136、OEM-B137



Outline dimensional drawing

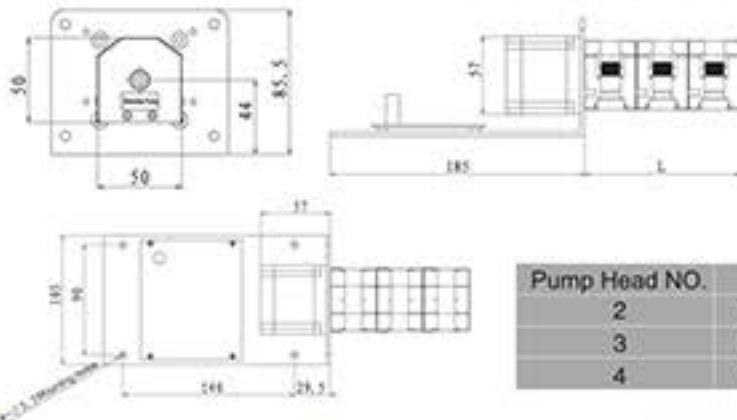

Pump Head NO.	L
2	75.6
3	113.4
4	151.2

OEM-B136 Technical Specification

Speed range	0.1–300 rpm
Flow rate range	0.0024–190mL/min
Speed control	Preset speed and manual control
Control mode	External passive signal level (normally closed or open) mode; control the start/stop of direction
Display (optional)	0.96" OLED screen
Motor type	57 stepper motor
Input voltage	DC 24V
Power	< 60W

OEM-B137 Technical Specification





Speed range	0.1–300 rpm
Flow rate range	0.0024–190mL/min
Speed control	Pulse frequency speed control
Logic input voltage	5V, 12V, 24V(default 5V)
Multiple segments	1/16(recommend), 1, 1/2, 1/4, 1/8, 1/32, 1/64, 1/128
Motor type	57 stepper motor
Input voltage	DC 24V
Power	< 60W

Motor speed(rpm)	Dual channel Flow rates (mL/min) and tubing (1 mm wall thickness)			
	ID 1mm	ID 2mm	ID 2.5mm	ID 3mm
5	0.25	0.90	1.28	1.78
10	0.49	1.84	2.59	3.52
15	0.74	2.76	3.85	5.30
25	1.25	4.57	6.43	8.80
40	2.00	7.31	10.25	14.02
50	2.52	9.14	12.78	17.47
75	3.76	13.64	19.04	26.12
100	5.03	18.17	25.42	35.34
200	10.12	36.45	51.27	71.44
250	12.63	45.75	64.34	88.35
300	15.01	54.63	76.84	108.39

Motor speed(rpm)	Single channel Flow rates (mL/min) and tubing (1.6 mm wall thickness)			
	13" ID 0.8mm	14" ID 1.6mm	16" ID 3.1mm	25" ID 4.8mm
5	0.12	0.56	1.97	3.26
10	0.26	1.15	3.90	6.41
15	0.39	1.74	5.80	9.52
25	0.66	2.88	9.59	15.74
40	1.08	4.59	15.29	25.15
50	1.34	5.69	18.96	31.14
75	2.00	8.54	28.36	46.55
100	2.73	11.37	37.66	61.98
200	5.46	22.76	76.02	124.14
250	6.81	28.40	95.08	153.63
300	8.28	33.88	114.31	190.00



OEM-B130 Patent No. ZL 2013 2 0435227.X

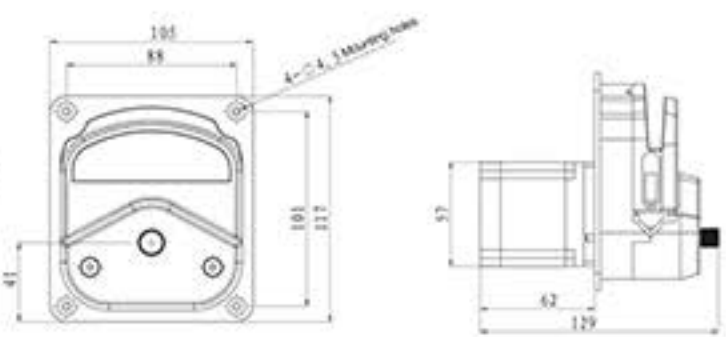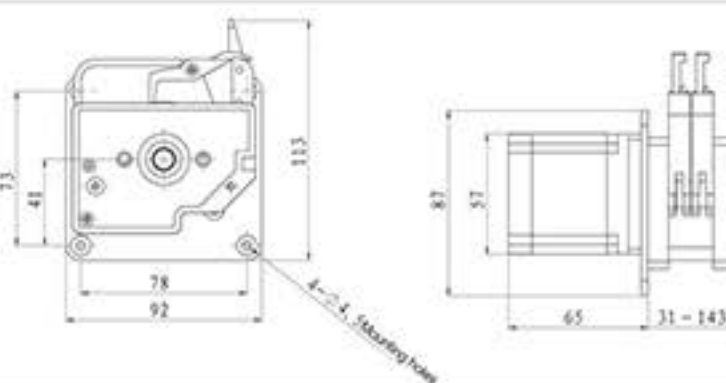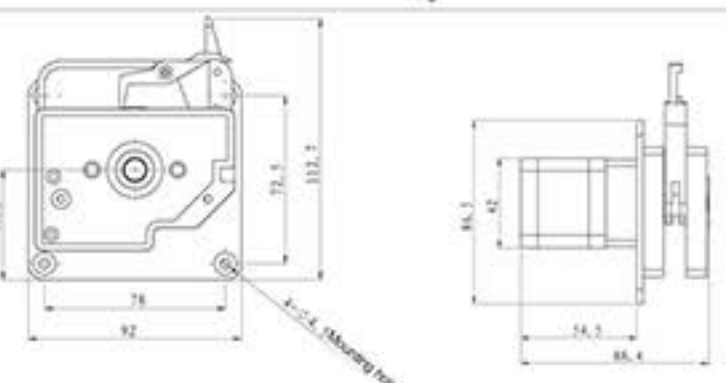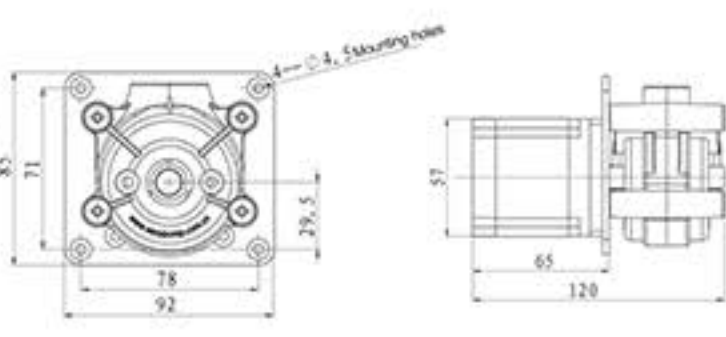 57 stepper motor	Outline dimensional drawing	
 57 stepper motor	Outline dimensional drawing	
 42 stepper motor	Outline dimensional drawing	
 57 stepper motor	Outline dimensional drawing	

Speed range	≤350 rpm
Max flow rate	1330mL/min
Motor type	57 stepper motor/42 stepper motor
Motor phase number	2
Motor step angle	1.8°
Motor phase voltage	3.75V
Motor phase current	3.0A

Tubing&Flow rate


Pump head	Tubing	Flow rate (mL/min)
YZ1515x	13', 14', 19', 16', 25', 17', 18'	≤1330
YZ2515x	15', 24'	≤1015
MC1-MC8	ID: 0.13-3.17mm Wall thickness: 0.8-1mm	≤49(10 roller)、≤65(6 roller)
SN Series	14', 16', 24'	≤875



OEM-B15 Series



Outline dimensional drawing

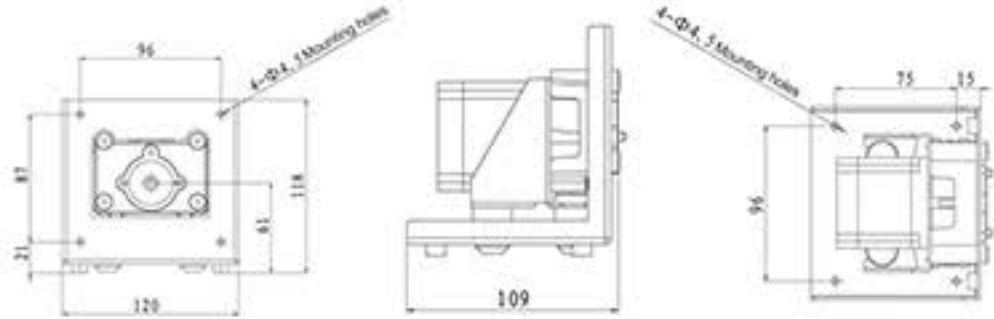

Suitable Pump Head



MC1~MC8



YZ1515x, YZ2515x



SN15, SN25

Speed range	≤350 rpm
Max flow rate	1330mL/min
Motor type	57 stepper motor
Motor phase number	2
Motor step angle	1.8°
Motor phase voltage	3.75V
Motor phase current	3.0A

Tubing&Flow rate

Pump head	Tubing	Flow rate (mL/min)
YZ1515x	13", 14", 19", 16", 25", 17", 18"	≤1330
YZ2515x	15", 24"	≤1015
MC1-MC8	ID: 0.13-3.17mm Wall thickness: 0.8-1mm	≤32(10 rollers), ≤45(6 rollers)
SN Series	14", 16", 24"	≤875

User driven / speed type control board



PD03 Main board technical specification

Speed range	0.1-350 rpm
Speed control	Pulse frequency speed control
Logic input voltage	5V, 12V, 24V(default 5V)
Multiple segments	1/16(recommend), 1, 1/2, 1/4, 1/8, 1/32, 1/64, 1/128
Input voltage	DC 24V
Power	<60W
Dimension (L*W)	98.48×82.68 mm



PMD07 Main board technical specification




Speed range	0.1-350 rpm
Speed control	Preset speed or manual control
Control mode	External passive signal level (normally closed or open) mode; control the start/stop and direction
Display (optional)	0.96" OLED screen
Input voltage	DC 24V
Power	<60W
Dimension (L*W)	98.48×82.68 mm



OEM Drive



Outline dimensional drawing

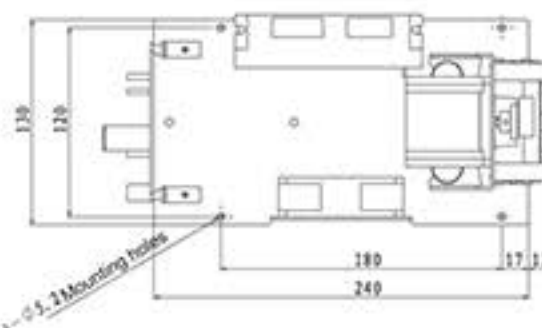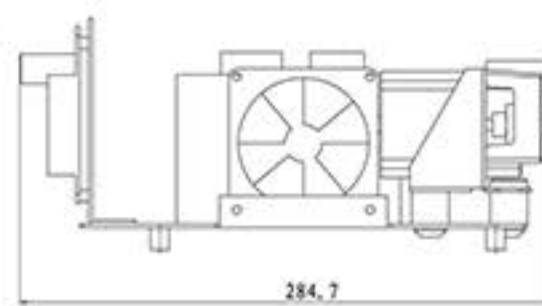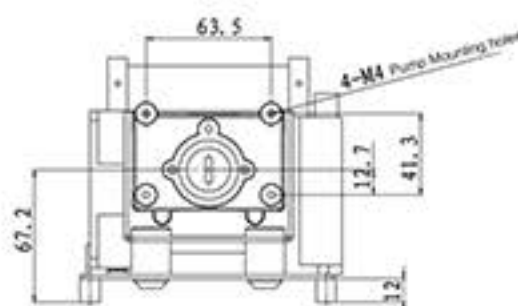

Suitable Pump Head



MC1~MC8



YZ1515x, YZ2515x



SN15, SN25

Pm01 Mainboard technical specification

Speed range	0.1~600 rpm(optional)
Flow rate range	0.000166~2280mL/min
Speed resolution	0.1rpm
Display	0.96" OLED screen(optional)
Speed control	Preset speed or membrane keypad speed control (optional)
Control mode	membrane keypad
External control	Start/stop,change direction signal ;Active switch signal: 5V (standard)
Communication interface	RS485 support Modbus protocol (RTU mode)
External control speed	0~5V、0~10V、4~20mA(optional)
Motor type	57 stepper motor
Input voltage	24V
Power	< 60W

Pm02 Mainboard technical specification

Speed range	0.1~600 rpm(optional)
Flow rate range	0.000067~2280mL/min
Speed resolution	0.1rpm
Display	0.96" OLED screen(optional)
Speed control	Preset speed or membrane keypad speed control (optional)
Control mode	Toggle switch
External control	Start/stop;Active switch signal: 5V (standard)
Communication interface	RS485 support Modbus protocol (RTU mode)
External control speed	0~5V、0~10V、4~20mA(optional)
Motor type	57 stepper motor
Input voltage	24V
Power	< 60W

Tubing&Flow rate

OEM Drive		Pump head & Flow rates (mL/min)						
		YZ1515x	YZ2515x	MC1~MC12(10)	MC1~MC12(6)	SN25	SN15-14	SN15-16
OEM Number (main board)&rotation speed	Tubing	13', 14', 19', 16' 25', 17', 18'	15', 24'	Inner diameter: 0.13~3.17mm Wall thickness: 0.86~1mm		24'	14'	16'
OEM-B145(PM01)	0.1~150rpm	0.007~570	0.17~435	0.000166~49 (suggest working rotating speed ≤150rpm)	0.000185~65 (suggest working rotating speed ≤150rpm)	0.25 ~ 375	0.024 ~ 36	0.088 ~ 132
OEM-B147(PM02)								
OEM-B146(PM01)	0.1~350rpm	0.007~1330	0.17~1015			0.25 ~ 875	0.024 ~ 84	0.088 ~ 308
OEM-B148(PM02)								
OEM-B221(PM01)	0.1~600rpm	0.007~2280	0.17~1740			0.25 ~ 1500	0.024 ~ 144	0.088 ~ 528
OEM-B222(PM02)								




Application case



OEM-Z Series



Outline dimensional drawing

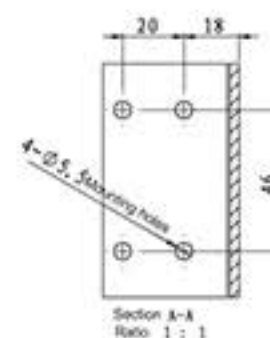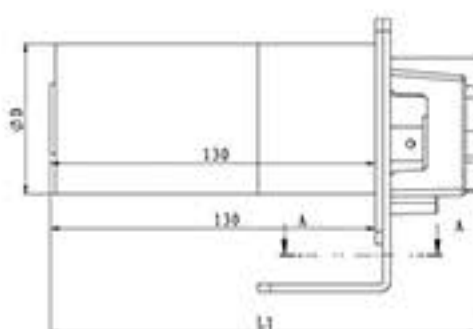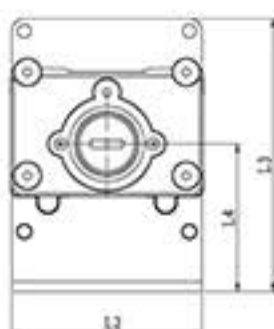

Suitable Pump Head



MC1~MC12



YZ1515x, YZ2515x



SN15, SN25

Product model	Install method	Speed (rpm)	Speed error	Suitable power	Motor type	Work environment		
OEM-Z101	L shape connection plate	60	± 10%	DC 12V/15W	Ø 60 DC geared motor	0-40℃		
OEM-Z102	L shape connection plate	100						
OEM-Z103	L shape connection plate	166						
OEM-Z104	L shape connection plate	200						
OEM-Z105	L shape connection plate	240						
OEM-Z106	L shape connection plate	333						
OEM-Z107	L shape connection plate	400						
OEM-Z108	L shape connection plate	500						
OEM-Z109	L shape connection plate	600						
OEM-Z110	—————	60						
OEM-Z111	L shape connection plate	60						
OEM-Z112	L shape connection plate	100						
OEM-Z113	L shape connection plate	166						
OEM-Z114	L shape connection plate	200						
OEM-Z115	L shape connection plate	240						
OEM-Z116	L shape connection plate	333						
OEM-Z117	L shape connection plate	400						
OEM-Z118	L shape connection plate	500						
OEM-Z119	L shape connection plate	600						
OEM-Z201	L shape connection plate	60	± 10%	DC 12V/25W	Ø 70 DC geared motor			
OEM-Z202	L shape connection plate	100						
OEM-Z203	L shape connection plate	150						
OEM-Z204	L shape connection plate	200						
OEM-Z205	L shape connection plate	300						
OEM-Z206	L shape connection plate	450						
OEM-Z207	L shape connection plate	500						
OEM-Z208	L shape connection plate	600						
OEM-Z211	L shape connection plate	60						
OEM-Z212	L shape connection plate	100						
OEM-Z213	L shape connection plate	150						
OEM-Z214	L shape connection plate	200						
OEM-Z215	L shape connection plate	300						
OEM-Z216	L shape connection plate	450						
OEM-Z217	L shape connection plate	500						
OEM-Z218	L shape connection plate	600						
OEM-Z301	L shape connection plate	600		± 10%			DC 24V/30W	Ø 60 DC motor
OEM-Z401	L shape connection plate	500		± 10%			DC 24V/40W	Ø 80 DC motor

Tubing&Flow rate



Pump head	Tubing	Flow rate (mL/min)
YZ1515x	13', 14', 19', 16', 25', 17', 18'	< 2280
YZ2515x	15', 24'	< 1740
MC	ID: 0.13–3.17mm Wall thickness: 0.86–1mm	< 45
SN Series	14', 16', 24'	< 1500



OEM-J Series



Outline dimensional drawing

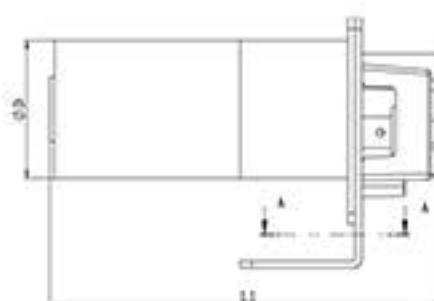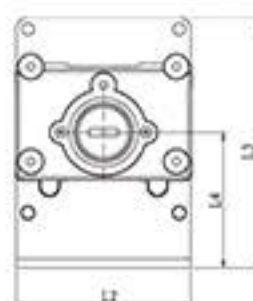

Suitable Pump Head



MC1~MC12



YZ1515x, YZ2515x



SN15, SN25



DZ25-3L



DZ25-6L

OEM-J

Product model	Install method	Speed (rpm)	Speed error	Suitable power	Motor type	Work environment
OEM-J01	L shape connection plate	60	± 10%	AC 220V/10W	Ø 60 AC geared motor	0-40℃
OEM-J02	L shape connection plate	100				
OEM-J03	L shape connection plate	166				
OEM-J04	L shape connection plate	200				
OEM-J05	L shape connection plate	300				
OEM-J06	L shape connection plate	417				
OEM-J08	L shape connection plate	60	± 10%	AC 220V/15W	Ø 70 AC geared motor	
OEM-J09	L shape connection plate	100				
OEM-J010	L shape connection plate	166				
OEM-J011	L shape connection plate	200				
OEM-J012	L shape connection plate	300				
OEM-J013	L shape connection plate	417				
OEM-J015	L shape connection plate	60	± 10%	AC 220V/25W	Ø 80 AC geared motor	
OEM-J016	L shape connection plate	120				
OEM-J017	L shape connection plate	300				
OEM-J018	L shape connection plate	500				
OEM-J019	L shape connection plate	60				
OEM-J020	L shape connection plate	120				
OEM-J021	L shape connection plate	300	± 10%	AC 220V/60W	Ø 90 AC geared motor	
OEM-J022	L shape connection plate	500				

OEM-JT

Product model	Install method	Speed range (rpm)	Speed error	Suitable power	Motor type	Work environment
OEM-JT01	L shape connection plate	0~60	± 10%	AC 220V/15W	Ø 70 AC geared motor	0~40℃
OEM-JT02	L shape connection plate	0~120				
OEM-JT03	L shape connection plate	0~300				
OEM-JT04	L shape connection plate	0~500				
OEM-JT05	L shape connection plate	0~60	± 10%	AC 220V/25W	Ø 80 AC geared motor	
OEM-JT06	L shape connection plate	0~120				
OEM-JT07	L shape connection plate	0~300				
OEM-JT08	L shape connection plate	0~500				
OEM-JT09	L shape connection plate	0~60	± 10%	AC 220V/60W	Ø 90 AC geared motor	
OEM-JT10	L shape connection plate	0~120				
OEM-JT11	L shape connection plate	0~300				
OEM-JT12	L shape connection plate	0~500				

Tubing&Flow rate

Pump head	Tubing	Flow rate (mL/min)
YZ1515x	13', 14', 19', 16', 25', 17', 18'	< 1900
YZ2515x	15', 24'	< 1450
MC	ID: 0.13~3.17mm Wall thickness: 0.8~1mm	< 45
SN Series	14', 16', 24'	< 1250
DZ25-3L	15', 24', 35', 36'	< 3000
DZ25-6L	15', 24', 35', 36'	< 5000



Silicone Tubing

Platinum-Cured Silicone Tubing

Slightly clarity, smooth surface, low protein binding levels, fewer potential leachable.

Ideal for pharmaceutical and biotechnology use, suitable temperature range -51~238 °C.

Micro Flow Rates Tubing

Tubing sizes	0.13×0.86	0.5×0.86	0.86×0.86	1.52×0.86	2.06×0.86	2.79×0.86	1×1	2×1	3×1	2.4×0.8
Tubing cross sections (1:1)										
Wall thickness (mm)	0.86						1.0		0.8	
Inside diameter (mm)	0.13	0.5	0.86	1.52	2.06	2.79	1.0	2.0	3.0	2.4
Maximum pressure (Mpa)	Continuous						0.1			
	Intermittent						0.1			

Basic Flow Rates Tubing

Tubing sizes	13'	14'	19'	16'	25'	17'	18'	15'	24'	35'	36'
Tubing cross sections (1:1)											
Wall thickness (mm)	1.6							2.4			
Inside diameter (mm)	0.8	1.6	2.4	3.1	4.8	6.4	7.9	4.8	6.4	7.9	9.6
Maximum pressure (Mpa)	Continuous				0.14	0.1	0.07	0.17		0.14	
	Intermittent				0.24	0.14	0.1	0.27		0.24	

Industrial Tubing

Tubing sizes	26'	73'	82'	86'	90'	88'	92'
Tubing cross sections (1:1)							
Wall thickness (mm)	3.3			6.3		4.8	
Inside diameter (mm)	6.4	9.6	12.7	9.5	19	12.7	25.4
Maximum pressure (Mpa)	Continuous			0.25			
	Intermittent			0.3			



SAINT-GOBAIN Tubing Tygon, PharMed BPT, Norprene etc

	A Tygon3350	B Tygon R-3603	C Norprene Chemical	D PharMed	E Norprene A-60-F
					
Formulation	Tygon3350	Tygon R-3603	Norprene Chemical	PharMed	Norprene A-60-F
Application	Pharmaceutical, cosmetic, medical and auto-analysis application.	General laboratory, food & beverage, biopharmaceutical, analytical instruments.	Excellent for chemical processing and general industrial applications. Food and beverage applications where extractables are a concern.	Cell and tissue culture work and pharmaceutical uses. Also good for light-sensitive samples.	Ideal for the food, dairy and beverage.
Advantages	Ultra-smooth; minimizes bacterial growth. Good for mild to medium concentration bases, salts and alcohols; odorless, tasteless, and nontoxic. Transparent.	Inexpensive tubing for general lab application. Nonaging, nonoxidizing. Clear for easy flow monitoring. Handles virtually all inorganic chemicals. Low gas permeability. Smooth bore; good for viscous fluids. High dielectric constant.	Norprene thermoplastic elastomer outer jacket with chemically inert Tygon® 2075 inner bore for excellent chemical resistance. Plasticizer-free to guard against extractables. Long flex life. Opaque beige.	Great for tissue and cell work-nontoxic and nonhemolytic; long service life minimizes risk of fluid exposure; reduces tubing costs and pump downtime. Opaque to UV and visible light to protect light-sensitive fluids. Heat sealable, bondable, and formable. Extremely low gas permeability.	Heat, ozone, and UV light resistant. Nonaging; nonoxidizing; superior acid and alkali resistance. Opaque beige.
Application Suitability	—————	ACIDS GOOD ALKALIES GOOD ORGANIC SOLVENTS NO PRESSURE GOOD VACUUM GOOD VISCOUS FLUIDS EXCELLENT STERILE FLUIDS GOOD	—————	ACIDS GOOD ALKALIES GOOD ORGANIC SOLVENTS NO PRESSURE GOOD VACUUM EXCELLENT VISCOUS FLUIDS GOOD STERILE FLUIDS EXCELLENT	—————
Physical characteristics	—————	Thermoplastic. PVC-based material with plasticizer. Firm (stiff) material. Transparent, clear.	—————	Thermoplastic elastomer. Polypropylene-based material with USP mineral oil. Excellent tensile strength. Firm (stiff) material. Opaque, beige.	—————
Temp. range	-75 to 450°F (-60~232°C)	-58 to 165°F (-50~74°C)	-76 to 165°F (-60~74°C)	-60 to 270°F (-59~135°C)	-60 to 275°F (-51~135°C)
Meets classifications	FDA 21 CFR 177.2600 USP Class VI EP 3.1.9. Exceeds 3A standards Manufactured according to GMP.	FDA 21 CFR 175.300	None.	None.	FDA 21 CFR 177.2600 NSF listed (Standard 51) Manufactured according to GMP.
Cleaning/ Sterilization	Ethylene oxide gamma irradiation, or autoclave for 30 min, 15psi (1 bar).	Unaffected by commercial sanitizers (with recommended procedures). Sterilize with ethylene oxide (ETO) or autoclave. To autoclave: Coil loosely in nonlinting cloth or paper, autoclave at 121°C (250°F). 1KG/cm ² (15psi) for 30 minutes (tubing will appear milky); air dry at max 66°C (150°F) for 2 to 2 ½ hours until clear.	Sterilize with ethylene-oxide(ETO), autoclave or gamma irradiation up to 2.5Mrad. Repeated autoclaving will not affect overall life.	Autoclave, ethylene oxide, or gamma irradiation.	Autoclave.



SHENCHEN

Peristaltic Pump Tubing

	 F Norprene A-60-G	 G Tygon F-4040-A	 H Tygon LFL	 I TYGON 2475	 K Viton
Formulation	Norprene A-60-G	Tygon F-4040-A	Tygon LFL	TYGON 2475	Viton
Application	For applications requiring excellent chemical, heat, ozone, and ultra-violet (UV) light resistance.	Fuels and industrial lubricants-gasoline, kerosene, heating oils, cutting compounds, and glycol-based coolants. Resists most hydrocarbons.	General laboratory use, provides longer life with peristaltic tubing pumps.	Sensitive fluid transfer applications requiring high purity.	Acid and solvent transfer, high-temperature.
Advantages	Best choice for vacuum/pressure applications. Offers longest life with good flow consistency. Heat and ambient ozone resistant. Good resistance to acids/alkalies. Black color hides dirt and dust. Heat sealable, nonaging, and nonoxidizing. High dielectric constant.	Resists embrittlement and swelling, ozone-and UV-resistant, with low-extractability. Translucent yellow.	Longest life of all Tygon® peristaltic tubing (1000hrs). Nonaging, nonoxidizing. Clear for easy flow monitoring. Broad chemical resistance; low gas permeability. Smooth bore. Good for viscous fluids. High dielectric constant.	Plasticizer free, smooth inner surface (inhibits particulate buildup and bacterial growth), safely disposed of through incineration and nontoxic. Transparent.	The most chemical resistant tubing. Resistant to corrosives, solvents, and oils at elevated temperatures. Low gas permeability.
Application Suitability	ACIDS GOOD ALKALIES GOOD ORGANIC SOLVENTS NO PRESSURE EXCELLENT VACUUM EXCELLENT VISCOUS FLUIDS EXCELLENT STERILE FLUIDS NO	—————	ACIDS GOOD ALKALIES GOOD ORGANIC SOLVENTS NO PRESSURE GOOD VACUUM GOOD VISCOUS FLUIDS EXCELLENT STERILE FLUIDS POOR	—————	ACIDS EXCELLENT ALKALIES EXCELLENT ORGANIC SOLVENTS EXCELLENT PRESSURE GOOD VACUUM GOOD VISCOUS FLUIDS GOOD STERILE FLUIDS FAIR
Physical characteristics	Thermoplastic elastomer. Polypropylene-based material with USP mineral oil. Excellent tensile strength. Firm (stiff) material. Opaque, black. Manufactured according to GMP.	—————	Thermoplastic. PVC-based material with plasticizer. Firm (stiff) material. Transparent, clear.	—————	Thermal set rubber. Viton B (67% fluorine) Firm (stiff) material. Opaque, black. Manufactured according to GMP.
Temp. range	-60 to 270°F (-59~135°C)	-35 to 165°F (-37~74°C)	-58 to 165°F (-50~74°C)	-94 to 125°F (-70~52°C)	-25 to 400°F (-32~205°C)
Meets classifications	None.	Meets NSF-51 and 3A sanitary standards.	USP Class VI, FDA 21 CFR 175.300	FDA 21 CFR 177.1520, USP 23 Class VI, Manufactured according to GMP.	None.
Cleaning/ Sterilization	Sterilize by autoclave only.	Not recommended.	Sterilize by ETO/autoclave. Coil loosely in nonlinting cloth or paper; autoclave at 250°F(121°C), 15 psi (1kg/cm²), 30 minutes (tubing will appear milky); air dry at max 150°F (66°C) for 2 to 2 ½ hrs until clear.	Ethylene oxide or gamma irradiation.	Sterilization is not recommended.



A Filling nozzle



- When transferring or dispensing liquid reduce splashing and dropping off maximumly, most out of the splash reduction.
- Dimensional tolerance $\pm 0.05\%$, assure reproducibility.



Name	Material	Tube
Filling nozzle	304 stainless steel	13#/14#, 19#, 16#, 25#/15#, 17#/24#

B Filling Countersunk



Used for the output of tube, preventing the tube floating or absorbing the wall of container.

Name	Material	Tube
Counter sunk	304/316 stainless steel	13#, 14#, 19#, 16#, 25#, 17#, 18#, 15#, 24#, 35#, 36#, 26#, 73#, 82#

C Foot pedal switch



Used for the drive without in hand, all pumps with 1.8m power line.

D Optical fiber sensors



Optical fiber sensors with the feature of anti-electromagnetic interference and atomic radiation interference, mechanical property with thin diameter, soft, light weight; electric property with insulation, noninductivity; chemical properties with water resistance, corrosion resistance, high temperature resistance, can be connected with filling system, as input signal of missing bottles.

E Benchtop tubing cutter

Stainless steel blade, makes right-angle cuts in several sizes of plastic tubing.



F Hall sensor



When applied in the dispensing line, it can detect whether there is filling bottle in the production line. When the bottle approaches the sensor side, the switch action will be made without any mechanical contact or pressure, thereby providing filling control order to the pump. In the same way, when no filling bottle is detected, the stop filling control order is provided to the pump.

F PLC



HMI true color touch screen. Our company supplies standard MODBUS communication protocol, can connect with various HMI products perfectly. It has a mature system for supporting industrial equipment, handles well, strong logical function, widely used for various industrial control occasions.

H Tube connector



1 Straight tube connector

Tube
2.5 x 1; 3 x 1; 16#



2 "Y" tube connector

Tube
2 x 1; 19#



3 Tee tube connector

Diameter		Tube
inch	mm	
1/16"	1.6	13' 14'
1/8"	3.2	16'
3/16"	4.8	15' 25'
1/4"	6.4	17' 24'
3/8"	9.6	18' 35' 36' 73' 86'
1/2"	12.7	82' 88'



4 Reducer tube connector

Diameter1		Diameter2		Tube1	Tube2
inch	mm	inch	mm		
1/16"	1.6	1/8"	3.2	13' 14'	16'
1/8"	3.2	3/16"	4.8	16'	15' 25'
1/8"	3.2	1/4"	6.4	16'	17' 24'
1/4"	6.4	3/8"	9.6	17' 24'	18' 35' 36' 73' 86'
3/8"	9.6	1/2"	12.7	18' 35' 36' 73' 86'	82' 88'
1/2"	12.7	3/4"	19	82' 88'	90'
1/2"	12.7	1"	25.4	82' 88'	92'

I Support stand





Gear Pump

LabGV, LabGF, LabGN, LabGM

NEW



Model Number

■ LabGN ■ LabGV ■ LabGF ■ LabGM

Typical Application

■ Mainly for laboratory.

Features

- With ABS engineering plastic housing, streamlined body design.
- Brushless DC motor with high efficiency and maintenance free features.
- The pump head is compact type, stainless steel material, driven by magnetic, low noise, no pulsation.
- Suitable for transferring viscous, high temperature and high pressure liquid.

Technical Specifications

Speed range	100~3000 rpm
Testing time range	0.1s~9999h
Outlet pressure	≤0.5Mpa
Motor type	DC brushless motor
Can transfer liquid which include granules	Granule ID ≤ 20 μm
Can transfer viscous liquid	Viscosity ≤ 200cSt
External speed control signal	0~5V, 0~10V, 4~20mA for option
Start/stop	Passive switch signal, such as foot pedal Active switch signal: 5V, 12V, 24V for option
Communication interface	RS232, RS485 support MODBUS protocol(RTU mode)
Power supply	AC 220V ± 10% 50Hz/60Hz(standard) AC 110V ± 10% 50Hz/60Hz(optional)
Drive dimension	261.4 × 157.3 × 236.9mm
Drive weight	4.40kg
Power consumption	<100W
Condition temperature	0~40℃
Relative humidity	<80%
IP rate	IP31

Technical Specifications

Drive	Model	Display	Operating mode	Function	Motor Status Output Interface
LabGV	Intelligent flow rate type	4.3"-industrial-true color LCD Screen	Touch screen + Membrane keypad	Flow rate/motor speed display, Fixed time and volume, Fixed volume measurement, Timer start and stop functions.	open-collector output
LabGF	Intelligent dispensing type	4.3"-industrial-true color LCD Screen	Touch screen + Membrane keypad	Input dispensing volume, dispensing time, copy numbers, for accuracy filling	open-collector output
LabGN	Standard type	3.2"LCD Screen	Membrane keypad	Motor speed and setting datas display in the same screen, used for common transferring, support simple dispensing(0.1s~9999h)	open-collector output
LabGM	Basic type	4 digital LED	Membrane keypad	Display motor speed, used for common transferring, support for simple dispensing (0.5s~999s)	No

Drive	Pump head	Gear material	Flow rates mL/min	Output rated pressure(Mpa)	Fluid temperature
LabGV LabGF LabGN LabGM	MG204XD0PT0000	PEEK	27.7 ~ 830	0.5	-30 ~ 100(℃)
	MG204XD0TT0000	PTFE	27.7 ~ 830		-30 ~ 50(℃)
	MG209XD0PT0000	PEEK	56.67 ~ 1700		-30 ~ 100(℃)
	MG209XD0TT0000	PTFE	56.67 ~ 1700		-30 ~ 50(℃)
	MG213XD0PT0000	PEEK	80 ~ 2400		-30 ~ 100(℃)
	MG213XD0TT0000	PTFE	80 ~ 2400		-30 ~ 50(℃)



Syringe Pump



SPLab02-Dual channel

Model Number

- SPLab01 (1 channel)
- SPLab02 (2 channels)
- SPLab04 (4 channels)
- SPLab06 (6 channels)
- SPLab08 (8 channels)
- SPLab10 (10 channels)
- SPLab12 (12 channels)



SPLab01-Single channel

Features

The syringe pump with infusion/withdrawal mode, stainless steel body, and large-screen chromatic LCD is easy to operate. It can install and fix different types of injectors or syringes in the same time, and it is suitable for transferring micro flow rates with high precision (non-medical).

Technical Specifications

Working mode(6 kinds)	Infusion, withdrawal, infusion/withdrawal, withdrawal/infusion, continuous, additional mode (infusion/withdrawal, then exclude bubble).						
Channel number	1	2	4	6	8	10	12
Syringe size	10 μ L-100mL		10 μ L-10mL				
Flow rates	0.001 μ L/min-127mL/min		0.001 μ L/min-21.997mL/min				
Infusion volume per microstep	0.078 μ m						
Min. linear rate	5 μ m/min						
Max. linear rate	132mm/min						
Min. step rate	0.035ms/Mirco step						
Max. step rate	0.937sec/Mirco step						
Max.linear force	10kgf						
Advance per microstep	0.078 μ m						
Accuracy	≤ ± 0.5%						
Operating mode	Rotary encoded switch and membrane keypad						
Display mode	320×240TFT-LCD						
Power-off memory	Display the previous data parameter after power supply again						
Power supply	AC 220V ± 10%(standard),AC 110V ± 10%(optional)						
Operating condition	Temperature 0-40℃. Relative humidity <80%						




SPLab10-10 channels

Syringe	Inside diameter(mm)	Flow rates(μ L/min–mL/min)	
		SPLab01/SPLab02	SPLab03-SPLab10
50 μ L	1.03	0.004~0.110	0.004~0.110
1mL	4.61	0.083~2.203	0.083~2.203
2.5mL	7.28	0.208~5.492	0.208~5.492
5mL	10.30	0.417~10.993	0.417~10.993
10 mL	14.57	0.834~21.997	0.834~21.997
25 mL	23.03	2.083~54.958	—
50 mL	32.57	4.166~109.921	—
60 mL	32.57	4.166~108.310	—


Installation Sizes Drawing (Unit: mm)

Channel No.	4	6	8	10	12
L	145	190	235	280	325


Syringe Pump External Control Schematic

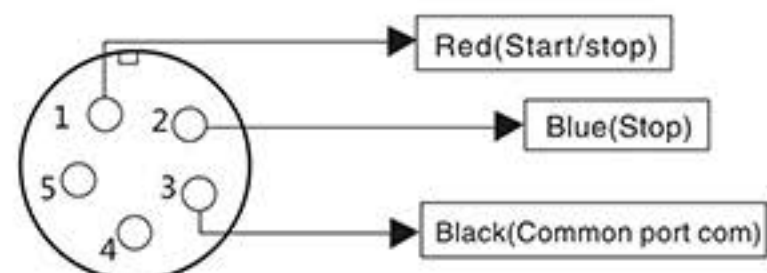


SPLab01

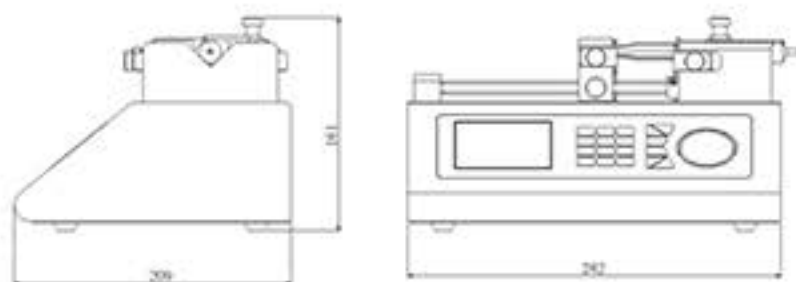


SPLab02

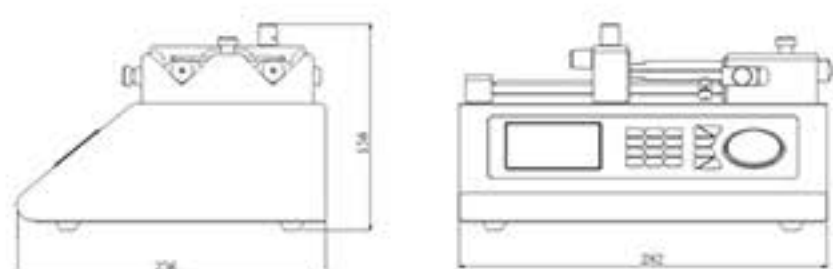

SPLab04

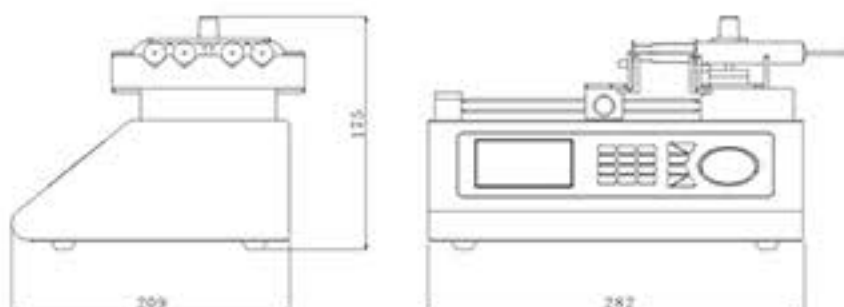

SPLab10

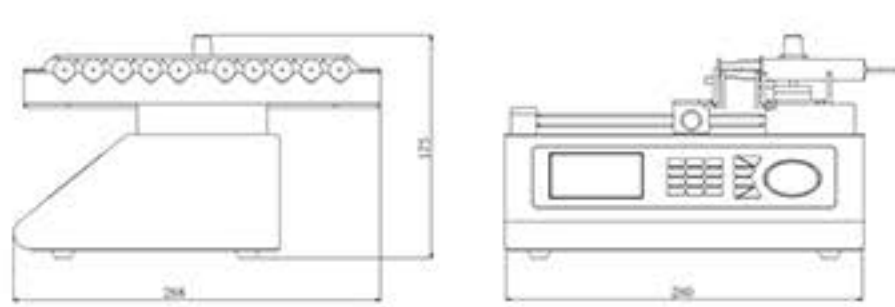



Shenchen precision pump
www.good-pump.com



Baoding Shencheng Precision Pump Co.,Ltd

Add: No.103,Building 2, ZhiDian Industrial Park,
FuXing East Road 999, Baoding,China.

Tel: 0086-312-6780680 6780681

Fax: 0086-312-6780636

Email: info@good-pump.com

Website: www.good-pump.com

